Database of rat liver proteins

N. Leigh Anderson Ricardo Esquer-Blasco Jean-Paul Hofmann Norman G. Anderson

Large Scale Biology Corporation, Rockville, MD

Lal et al., 09/002,485, filed December 31, 1997 (PF-0459)

Exhibit "L" attached to Declaration of John C. Rockett, Ph.D.

A two-dimensional gel database of rat liver proteins useful in gene regulation and drug effects studies

A standard two-dimensional (2-D) protein map of Fischer 344 rat liver (F344MST3) is presented, with a tabular listing of more than 1200 protein species. Sodium dodecyl sulfate (SDS) molecular mass and isoelectric point have been established, based on positions of numerous internal standards. This map has been used to connect and compare hundreds of 2-D gels of rat liver samples from a variety of studies, and forms the nucleus of an expanding database describing rat liver proteins and their regulation by various drugs and toxic agents. An example of such a study, involving regulation of cholesterol synthesis by cholesterol-lowering drugs and a high-cholesterol diet, is presented. Since the map has been obtained with a widely used and highly reproducible 2-D gel system (the Iso-Dalt* system), it can be directly related to an expanding body of work in other laboratories

Contents

1	Introduction	907
	Material and methods	908
_	2.1 Sample preparation	908
	2.2 Two-dimensional electrophoresis	909
	2.3 Staining	909
	2.4 Positional standardization	909
	2.5 Computer analysis	909
	2.6 Graphical data output	910
	2.7 Experiment LSBC04	910
3	Results and discussion	910
_	3.1 The rat liver protein 2-D map	910
	3.2 Carbamylated charge standards computed prs	
	and molecular mass standardization	911
	3.3 An example of rat liver gene regulation: Chol-	
	esterol metabolism	911
	3.3.1 MSN 413 (putative cytosolic HMG-CoA	
	synthase) and sets of spots regulated co-	
	ordinately or inversely	911
	3.3.2 MSN 235 and corregulated spots	912
	3.3.3 An example of an anti-synergistic effect	912
	3.3.4 Complexity of the cholesterol synthesis	
	pathway	912
4	Conclusions	912
5	References	912
6	Addendum 1: Figures 1–13	914
7	Addendum 2: Tables 1-4	923
	Table 1. Master table of proteins in rat liver data-	
	base	923
	Table 2. Table of some identified proteins	928
	Table 3. Computed pl's of two sets of carbamylated	
	protein standards: rabbit muscle CPK and	
	human Hb	929
	Table 4. Computed pl's of some known proteins re-	
	lated to measured CPK pl's	930

Correspondence: Dr. N. Leigh Anderson, Large Scale Biology Corporation, 9620 Medical Center Drive, Rockville, MD 20850, USA

Abbreviations: CBB, Coomassie Brilliant Blue; CPK, creatine phosphokinase; 2-D, two-dimensional; IEF, isoelectric focusing; MSN, master spot number; NP-40, Nonidet P-40, SDS, sodium dodecyl sulfate

1 Introduction

High-resolution two-dimensional electrophoresis of proteins, introduced in 1975 by O'Farrell and others [1-4], has been used over the ensuing 16 years to examine a wide variety of biological systems, the results appearing in more than 5000 published papers. With the advent of computerized systems for analyzing two-dimensional (2-D) gel images and constructing spot databases, it is also possible to plan and assemble integrated bodies of information describing the appearance and regulation of thousands of protein gene products [5, 6]. Creating such databases involves amassing and organizing quantitative data from thousands of 2-D gels, and requires a substantial commitment in technology and resources.

Given the long-term effort required to develop a protein database, the choice of a biological system takes on considerable importance. While in vitro systems are ideal for answering many experimental questions, especially in cancer research and genetics, our experience with cell cultures and tissue samples suggests that some in vivo approaches could have major advantages. In particular, we have noticed that liver tissue samples from rats and mice appear to show greater quantitative reproducibility (in terms of individual protein expression) than replicate cell cultures. This is perhaps a natural result of the homeostasis maintained in a complete animal vs. the well-known variability of cell cultures, the latter due principally to differences in reagents (e.g., fetal bovine serum), conditions (e.g., pH) and genetic "evolution" of cell lines while in culture. It is also more difficult to generate adequate amounts of protein from cell culture systems (particularly with attached cells), forcing the investigator to resort to radioisotope-based or silver-based staindetection methods. While these methods are more sensitive (sometimes much more sensitive) than the Coomassie Brilliant Blue (CBB) stain typically used for protein detection in "large" protein samples, they are generally more variable, more labor-intensive and, in the case of radiographic methods, may generate highly "noisy" images, due to the properties of the films used. By contrast, large protein samples can easily be prepared from liver using urea/Nonidet P-40 (NP-40) solubilization and stained with CBB, which has the advantage of being easily reproducible [8]. Finally, there remains the question of the "truthfulness" of many in vitro systems as compared to their in vivo analogs; how great are the changes caused by the introduction into a culture and the associated shift to strong selection for growth, and how do these affect experimental outcomes? Hence the apparent advantages of in vitro systems, in terms of experimental manipulation, may be counterbalanced by other factors relating to 2-D data quality.

There is a second important class of reasons for exploring the use of an in vivo biological system such as the liver. Historically, there have been two broad approaches to the mechanistic dissection of biochemical processes in intact cellular systems: genetics (a search for informative mutants) and the use of chemical agents (drugs and chemical toxins). Both approaches help us to understand complex systems by disrupting some specific functional element and showing us the result. With the development of techniques for genetic manipulation and cloning, the genetic approach can be effectively applied either in vitro or in vivo, although the in vitro route is usually quicker. The chemical approach can also be applied to either sort of biological system; here, however, the bulk of consistently acquired information is in experimental animals (rats and mice). While most biologists know a short list of compounds having specific, experimentally useful effects (e.g., inhibitors of protein synthesis, ionophores, polymerase inhibitors, channel blockers, nucleotide analogs, and compounds affecting polymerization of cytoskeletal proteins), there is a much larger number of interesting chemically-induced effects, most of them characterized by toxicologists and pharmacologists in rodent systems. Just as a thorough genetic analysis would involve saturating a genome with mutations, it is possible to imagine a saturating number of drugs, the analysis of whose actions would reveal the complete biochemistry of the cell. While organized drug discovery efforts usually target specific desired effects, the nature of the process, with its dependence on screening large numbers of compounds, necessarily produces many unanticipated effects. It is therefore reasonable to suppose that the required broad range of compounds necessary to achieve "biochemical saturation" may be forthcoming; in fact, it may already exist among the hundreds of thousands of compounds that failed to qualify as drugs.

Among organs, the liver is an obvious choice for the study of chemical effects because of its well-known plasticity and responsiveness. The brain appears to be quite plastic (e.g. [7]), but it is a complicated mixture of cell types requiring skillful dissection for most experiments. The kidney, while quite responsive, also presents a potentially confounding mixture of cell types. The liver, by contrast, is made up of one predominant cell type which is easy to solubilize: the hepatocyte, representing more than 95% of its mass. Most importantly, the liver performs many homeostatic functions that require rapid modulation of gene expression. It appears that most chemical agents tested affect gene expression in the liver at some dosage (N. Leigh Anderson, unpublished observations), an interesting contrast to our earlier work with lymphocytes, for example, which seem to be much less responsive. Such results conform to the expectation that cells with a homeostatic, physiological role should be more plastic than cells differentiated for a purpose dependent on the action of a limited number of specific genes.

The liver also allows the parallels between in vitro and in vivo systems to be examined in detail. Significant progress

has been made in the development of mouse, rat and human hepatocyte culture systems, as well as in precision-cut tissue slices. Using such an array of techniques, it is possible to assemble a matrix of mammalian systems including mouse and rat in vivo on one level and mouse, rat and human in vitro on a second level, and to compare effects between species and between systems. This approach allows us to draw informed conclusions regarding the biochemical "universality" of biological responses among the mammals, and to offer some insight into the validity of in vitro approaches for toxicological screening. We believe this data will be necessary if in vitro alternatives are to achieve wide usage in government-mandated safety testing of drugs, consumer products and industrial and agricultural chemicals.

A number of interesting studies have been published using 2-D mapping to examine effects in the rodent liver. A number of investigarors have made use of the technique to screen for existing genetic variants [8–11] or induced mutations [12–14], mainly in the mouse. This work builds on the wealth of genetic information available on the mouse and its established position as a mammalian mutation-detection system. While some studies of chemical effects have been undertaken in the mouse [15–17], most have used the rat [18–23]. The examination of the cytochrome p-450 system, in particular, has been carried out almost exclusively on the rat [24, 25].

These considerations lead us to conclude that rodent liver offers the best opportunity to systematically examine an array of gene regulation systems, and ultimately to build a predictive model of large-scale mammalian gene control. The basic underlying foundation of such a project is a reliable, reproducible master 2-D pattern of liver, to which ongoing experimental results can be referred. In this paper, we report such a master pattern for the acidic and neutral proteins of rat liver (pattern F344MST3). In future, this master will be supplemented by maps of basic proteins, and analogous maps of mouse and human liver.

2 Materials and methods

2.1 Sample preparation

Liver is an ideal sample material for most biochemical studies, including 2-D analysis. A sample is taken of approximately 0.5 g of tissue from the apical end of the left lobe of the liver. Solubilization is effected as rapidly as practical; a delay of 5–15 min appears to cause no major alteration in liver protein composition if the liver pieces are kept cold (e.g., on ice) in the interim. In the solubilization process, the liver sample is weighed, placed in a glass homogenizer (e.g., 15 mL Wheaton); 8 volumes of solubilizing solution*

^{*} The solubilizing solution is composed of 2% NP-40 (Sigma), 9 M urea (analytical grade, e.g., BDH or Bio-Rad), 0.5% dithiothreitol (DTT; Sigma) and 2% carrier ampholytes (pH 9-11 LKB: these come as a 20% stock solution, so 2% final concentration is achieved by making the final solution 10% 9-11 Ampholine by volume). A large batch of solubilizer (several hundred mL) is made and stored frozen at -80°C in aliquots sufficient to provide enough for one day's estimated sample preparation requirement. The solution is never allowed to become warmer than room temperature at any stage during preparation or thawing for use, since heating of concentrated urea solutions can produce contaminants that covalently modify proteins producing artifactual charge shifts. Once thawed, any unused solubilizer is discarded.

is added (i.e., 4 mL per 0.5 g tissue) and the mixture is homogenized using first the loose- and then then the tight-fitting glass pestle. This takes approximately 5 strokes with each pestle and is carried out at room temperature because urea would crystallize out in the cold. Once the liver sample is thoroughly homogenized in the solubilizer, it is assumed that all the proteins are denatured (by the chaotropic effect of the urea and NP-40 detergent) and the enzymes inactivated by the high pH (~9.5). Therefore these samples may be kept at room temperature until they can be centrifuged or frozen as a group (within several hours of preparation). The samples are centrifuged for 6×10^6 g min (e.g., 500 000 × g for 12 min using a Beckman TL-100 centrifuge). The centrifuge rotor is maintained at just below room temperature (e.g., 15-20°C), but not too cold, so as to prevent the precipitation of urea. The centrifuge of choice is a Beckman TL-100 because of the sample tube sizes available, but any ultracentrifuge accepting smallish tubes will suffice. When an appropriate centrifuge is not available near the site of sample preparation, samples can be frozen at -80°C and thawed prior to centrifugation and collection of supernatants. Each supernatant is carefully removed following centrifugation and aliquoted into at least 4 clean tubes for storage. This is done by transferring all the supernatant to one clean tube, mixing this gently (to assure homogeneous composition) and then dividing it into 4 aliquots. The aliquots are frozen immediately at -80°C. These multiple aliquots can provide insurance against a failed run or a freezer breakdown.

2.2 Two-dimensional electrophoresis

Sample proteins are resolved by 2-D electrophoresis using the 20 × 25 cm Iso-Dalt[®] 2-D gel system ([26-29]; produced by LSB and by Hoefer Scientific Instruments, San Francisco) operating with 20 gels per batch. All first-dimensional isoelectric focusing (IEF) gels are prepared using the same single standardized batch of carrier ampholytes (BDH 4-8A in the present case, selected by LSB's batchtesting program for rat and mouse database work**). A 10 μL sample of solubilized liver protein is applied to each gel, and the gels are run for 33 000 to 34500 volt-hours using a progressively increasing voltage protocol implemented by a programmable high-voltage power supply. An Angelique™ computer-controlled gradient-casting system (produced by LSB) is used to prepare second-dimensional sodium dodecyl sulfate (SDS) polyacrylamide gradient slab gels in which the top 5% of the gel is 11%T acrylamide, and the lower 95% of the gel varies linearly from 11% to 18%T.

This system has recently been modified so as to employ a commercially available 30.8%T acrylamide/N,N-methylenebisacrylamide prepared solution (thus avoiding the handling of the solid acrylamide monomer) and three additional stock solutions: buffer (made from Sigma pre-set Tris), persulfate and N,N,N',N-tetramethylethylenediamine (TEMED). Each gel is identified by a computer-printed filter paper label polymerized into the lower left corner of the gel. First-dimensional IEF tube gels are loaded

directly (as extruded) onto the slab gels without equilibration, and held in place by polyester fabric wedges (Wedgies^m, produced by LSB) to avoid the use of hot agarose. Second-dimensional slab gels are run overnight, in groups of 20, in cooled DALT tanks (10°C) with buffer circulation. All run parameters, reagent source and lot information, and notations of deviation from expected results are entered by the technician responsible on a detailed, multi-page record of the experiment.

2.3 Staining

Following SDS-electrophoresis, slab gels are stained for protein using a colloidal Coomassie Blue G-250 procedure in covered plastic boxes, with 10 gels (totalling approximately 1 L of gel) per box. This procedure (based on the work of Neuhoff [30, 31]) involves fixation in 1.5 L of 50% ethanol and 2% phosphoric acid for 2 h, three 30 min washes, each in 2L of cold tap water, and transfer to 1.5L of 34% methanol, 17% ammonium sulfate and 2% phosphoric acid for 1 h, followed by the addition of a gram of powdered Coomassie Blue G-250 stain. Staining requires approximately 4 days to reach equilibrium intensity, whereupon gels are transferred to cool tap water and their surfaces rinsed to remove any particulate stain prior to scanning. Gels may be kept for several months in water with added sodium azide. The water washes remove ethanol that would dissolve the stain (and render the system noncolloidal, with high backgrounds). The concentrated ammonium sulfate and methanol solution is diluted by equilibration with the water volume of the gels to automatically achieve the correct final concentrations for colloidal staining. Practical advantages of this staining approach can be summarized as follows: (i) the low, flat background makes computer evaluation of small spots (max OD < 0.02) possible, especially when using laser densitometry; (ii) up to 1500 spots can be reliably detected on many gels (e.g., rat liver) at loadings low enough to preserve excellent resolution; and (iii) reproducibility appears to be very good: at least several hundred spots have coefficients of reproducibility less than 15%. This value is at least as good as previous CBB methods, and significantly better than many silver stain systems.

2.4 Positional standardization

The carbamylated rabbit muscle creatine phosphokinase (CPK) standards [32] are purchased from Pharmacia and BDH. Amino acid compositions, and numbers of residues present in proteins used for internal standardization, are taken from the Protein Identification Resource (PIR) sequence database [33].

2.5 Computer analysis

Stained slab gels are digitized in red light at 134 micron resolution, using either a Molecular Dynamics laser scanner (with pixel sampling) or an Eikonix 78/99 CCD scanner. Raw digitized gel images are archived on high-density DAT tape (or equivalent storage media) and a greyscale videoprint prepared from the raw digital image as hard-copy backup of the gel image. Gels are processed using the Kepler® software system (produced by LSB), a commercially available workstation-based software package built on

This material (succeeding certified batches of which are available from Hoefer Scientific Instruments) has the most linear pH gradient produced by any ampholyte tested except for the Pharmacia wide range (which has an unacceptable tendency to bind high-molecular weight acidic proteins, causing them to streak).

some of the principles of the earlier TYCHO system [34–41]. Procedure PROC008 is used to yield a spotlist giving position, shape and density information for each detected spot. This procedure makes use of digital filtering, mathematical morphology techniques and digital masking to remove the background, and uses full 2-D least-squares optimization to refine the parameters of a 2-D Gaussian shape for each spot. Processing parameters and file locations are stored in a relational database, while various log files detailing operation of the automatic analysis software are archived with the reduced data. The computed resolution and level of Gaussian convergence of each gel are inspected and archived for quality control purposes.

Experiment packages are constructed using the Kepler experiment definition database to assemble groups of 2-D patterns corresponding to the experimental groups (e.g., treated and control animals). Each 2-D pattern is matched to the appropriate "master" 2-D pattern (pattern F344MST3 in the case of Fischer 344 rat liver), thereby providing linkage to the existing rodent protein 2-D databases. The software allows experiments containing hundreds of gels to be constructed and analyzed as a unit, with up to 100 gels displayed on the screen at one time for comparative purposes and multiple pages to accommodate experiments of > 1000 gels. For each treatment, proteins showing significant quantitative differences vs. appropriate controls are selected using group-wise statistical parameters (e.g., Student's t-test, Kepler® procedure STUDENT). Proteins satisfying various quantitative criteria (such as P< 0.001 difference from appropriate controls) are represented as highlighted spots onscreen or on computer-plotted protein maps and stored as spot populations (i.e., logical vectors) in a liver protein database. Quantitative data (spot parameters, statistical or other computed values) are stored as real-valued vectors in the database. Analysis of coregulation is performed using a Pierson product-moment correlation (Kepler procedure CORREL) to determine whether groups of proteins are coordinately regulated by any of the treatments. Such groups can be presented graphically on a protein map, and reported together with the statistical criteria used to assess the level of coregulation. Multivariate statistical analysis (e.g., principal components' analysis) is performed on data exported to SAS (SAS Institute).

2.6 Graphical data output

Graphical results are prepared in GKS and translated within Kepler® into output for any of a variety of devices. Linedrawing output is typically prepared as Postscript and printed on an Apple Laserwriter. Detailed maps presented here have been generated using an ultra-high-resolution Postscript-compatible Linotronic output device. Greyscale graphics are reproduced from the workstation screen using a Seikosha videoprinter. Patterns are shown in the standard orientation, with high molecular mass at the top and acidic proteins to the left.

2.7 Experiment LSBC04

In the study described here 12-week-old Charles River male F344 rats were used. Diets were prepared at LSB, based on a Purina 5755M Basal Purified Diet. Lovastatin and cholestyramine were obtained as prescription pharma-

ceuticals, ground and mixed with the diet at concentrations of 0.075% and 1%, respectively. The high cholesterol diet was Purina 5801M-A (5% cholesterol plus 1% sodium cholate in the control diet). Animal work was carried out by Microbiological Associates (Bethesda, MD). Animals were acclimatized for one week on the control diet, fed test or control diets for one week, and sacrificed on day 8. Average daily doses of lovastatin and cholestyramine in appropriate groups were 37 mg/kg/day and 5 g/kg/day, respectively, based on the weight of the food consumed. Liver samples were collected and prepared for 2-D electrophoresis according to the standard liver protocol (homogenization in 8 volumes of 9 m urea, 2% NP-40, 0.5% dithiothreitol, 2% LKB pH 9-11 carrier ampholytes, followed by centrifugation for 30 min at $80000 \times g$). Kidney, brain and plasma samples were frozen. Gels were run as described above, and the data was analyzed using the Kepler® system. Gels were scaled, to remove the effect of differences in protein loading, by setting the summed abundances of a large number of matched spots equal for each gel (linear scaling).

3 Results and discussion

3.1 The rat liver protein 2-D map

F344MST3 is a standard 2-D pattern of rat liver proteins, based on the Fischer 344 strain. This pattern was initiated from a single 2-D gel and extensively edited in an experiment comparing it to a range of protein loads, so as to include both small spots and well-resolved representations of high-abundance spots. More than 700 rat liver 2-D patterns have been matched to F344MST3 in a series of drug effects and protein characterization experiments, and numerous new spots (induced by specific drugs, for instance) have been added as a result. A modified version including additional spots present in the Sprague-Dawley outbred rat has also been developed (data not shown). Figure 1 shows a greyscale representation and Fig. 2 a schematic plot of the master pattern. More than 1200 spots are included, most of which are visible on typical gels loaded with 10 µL of solubilized liver protein prepared by the standard method and stained with colloidal Coomassie Blue. Master spot numbers (MSN's) have been assigned to all proteins, and appear in the following figures, each showing one quadrant of the pattern. Figure 3 shows the upper left (acidic, high molecular mass) quadrant, Fig. 4 the upper right (basic, high molecular mass) quadrant, Fig. 5 the lower left (acidic, low molecular mass) quadrant, and Fig. 6 the lower right (basic, low molecular mass) quadrant. The quadrants overlap as an aid to moving between them. The gel position (in 100 micron units), isoelectric point (relative to the CPK internal pI standards) and SDS molecular mass (from the calibration curve in Fig. 8) are listed for each spot (Table 1). Because of the precision of the CPK-pI values, these parameters can be used to relate spot locations between gel systems more reliably than using pI measurements expressed as pH. A major objective of current studies is the identification of all major spots corresponding to known liver proteins, as well as rigorous definitions of subcellular organelle contents. Of particular interest to us is the parallel development of identifications in the rat and mouse liver maps, allowing detailed comparisons of gene expression effects in the two systems. The results of these studies will be presented systematically in a later edition of this database, but we include here a useful series of 22 orienting identifications as an aid to other users of the rat liver pattern (Table 2).

3.2 Carbamylated charge standards, computed pPs and molecular mass standardization

We have previously shown that the use of a system of closely-spaced internal pI markers (made by carbamylating a basic protein) offers an accurate and workable solution to the problem of assigning positions in the pI dimension [32]. The same system, based on 36 protein species made by carbamylating rabbit muscle CPK, has been used here to assign pI's to most rat liver acidic and neutral proteins. The standards were coelectrophoresed with total liver proteins, and the standard spots added to a special version of the master pattern F344MST3. The gel X-coordinates of all liver protein spots lying within the CPK charge train were then transformed into CPK pI positions by interpolation between the positions of immediately adjacent standards (Table 1) using a Kepler® vector procedure.

It has proven possible to compute fairly accurate pI values for many proteins from the amino acid composition [42]. We have attempted here to test a further elaboration of this approach, in which we computed pl's for the CPK standards themselves, based on our knowledge of the rabbit muscle CPK sequence and the fact that adjacent members of the charge train typically differ by blockage of one additional lysine residue (Table 3). We compared these values to similar computed pl's for an additional set of carbamylated standards made from human hemoglobin beta chains and a series of rat liver and human plasma proteins of known position and sequence (Fig. 7, Table 4). The result demonstrates good concordance between these systems. Two proteins show significant deviations: liver fatty-acid binding protein (FABP; #1 in Table 4) and protein disulphide isomerase (#20 in the table). The FABP spot present on F344MST3 may represent a charge-modified version of a more basic parent spot closer to the expected pl, not resolved in the IEF/SDS gel. Of particular importance is the fact that, by comparing computed pl's of sequenced but unlocated proteins with the CPK pl's, we can assign a probable gel location without making any assumptions regarding the actual gel pH gradient. This offers a useful shortcut, given the vagaries of pH measurement on small diameter IEF gels. We have used this approach to compute the CPK pr of all rat and mouse proteins in the PIR sequence database, as an aid to protein identification (data not shown).

In order to standardize SDS molecular weight (SDS-MW), we have used a standard curve fitted to a series of identified proteins (Fig. 8). Rather than using molecular mass per se, we have elected to use the number of amino acids in the polypeptide chain, as perhaps a better indication of the length of the SDS-coated rod that is sieved by the second dimension slab. The resulting values were multiplied by 112 (the weighted average mass of amino acids in sequenced proteins) to give predicted molecular masses. Because we use gradient slabs, we have not constrained the fitted curve to conform to any predetermined model; rather we tried many equations and selected the best using the program "Tablecurve" on a PC. The equation chosen was y = a + bx + c/x, where y is the number of residues, x is the gel

Y coordinate, a is 511.83, b is -0.2731 and c is 33183801. The resulting fit appears to be fairly good over a broad range of molecular mass.

3.3 An example of rat liver gene regulation: Cholesterol metabolism

Experiment LSBC04 was designed as a small-scale test of the regulation of cholesterol metabolism in vivo by three agents included in the diet: lovastatin (Mevacor®, an inhibitor of HMG-CoA reductase); cholestyramine (a bile acid sequestrant that has the effect of removing cholesterol from the gut-liver recirculation); and cholesterol itself. The first two agents should lower available cholesterol and the third should raise it, allowing manipulation of relevant gene expression control systems in both directions. Such an experiment offers an interesting test of the 2-D mapping system since most of the pathway enzymes are present in low abundance, many are membrane-bound and difficult to solubilize, and the pathway itself is complex. Approximately 1000 proteins were separated and detected in liver homogenates. Twenty-one proteins were found to be affected by at least one treatment, and these could be divided into several coregulated groups.

3.3.1 MSN 413 (putative cytosolic HMG-CoA synthase) and sets of spots regulated coordinately or inversely

One group of spots (including a spot assigned to the cytosolic HMG-CoA synthase, MSN 413) showed the expected increase in abundance with lovastatin or cholestyramine, the synergistic further increase with lovastatin and cholestyramine, and a dramatic decrease with the high cholesterol diet. Spot number 413 is the most strongly regulated protein in the present experiment, showing a 5- to 10-fold induction after a 1 week treatment with 0.075 % lovastatin and 1% cholestyramine in the diet (Figs. 9 and 10). Its expression follows precisely the expectation for an enzyme whose abundance is controlled by the cholesterol level; it is progressively increased from the control levels by cholestyramine, lovastatin and lovastatin plus cholestyramine, and it sinks below the threshold of detection in animals fed the high cholesterol diet. This spot has been tentatively identified as the cytosolic HMG-CoA synthase, based on a reaction with an antiserum to that protein provided by Dr. Michael Greenspan at Merck Sharp & Dohme Research Laboratories. This enzyme lies immediately before HMG-CoA reductase in the liver cholesterol biosynthesis pathway, and is known to be co-regulated with it. Spot 413 has an SDS molecular weight of about 54 000 and a CPK pl of -11.4, in reasonably close agreement with a molecular weight of 57300 and a CPK pl of -15.7 computed from the known sequence of the hamster enzyme [43].

Using a classical product-moment correlation test (Kepler procedure CORREL), a series of five additional spots was found to be coregulated with 413. The level of correlation was exceedingly high (> 95%). Two of these, 1250 and 933, are at similar molecular weights and approximately one charge more acidic than 413 (Fig. 9), indicating that they may be covalently modified forms of the 413 polypeptide. This suspicion is strengthened by the observation that both spots are also stained by the antibody to cytosolic HMG-CoA synthase. The remaining three correlated spots appear

to comprise an additional related pair (1253 and 1001) of around 40 kDa and a single spot (1119) of around 28 kDa. Because these two presumed proteins are present at substantially lower abundances than 413, and because the cytosolic HMG-CoA synthase is reported to consist of only one type of polypeptide, they are likely to represent other, very tightly coregulated enzymes. A second group of six spots was selected based on a regulatory pattern close to the inverse of that for spot 413 (MSN's 34, 79, 178, 182, 204, 347; data not shown). For these proteins, the lowest level of expression occurs with exposure to lovastatin plus cholestyramine and the highest level upon exposure to the high-cholesterol diet. Spots 182 and 79 are highly correlated and lie about one charge apart at the same molecular weight; they may thus be isoforms of a single protein. The other four spots probably represent additional enzymes or subunits.

3.3.2 MSN 235 and coregulated spots

A third group of five spots, mainly comprised of mitochondrial proteins including putative mitochondrial HMG-CoA synthase spots, showed a modest induction by lovastatin alone, but little or no effect with any of the other treatments (including the combination of lovastatin and cholestyramine; Fig. 12). This result is intriguing because lovastatin was expected to affect only the regulation of enzymes of cholesterol synthesis, which is entirely extra-mitochondrial. Three of the spots (235, 134, 144) form a closelypacked triad at approximately 30 kDa, and are likely to represent isoforms of one protein. All three spots are stained by an antibody to the mitochondrial form of HMG-CoA synthase obtained from Dr. Greenspan. Subcellular fractionation indicates a mitochondrial location. The other two spots (633 at about 38 kDa and 724 at about 69 kDa) are each present at lower abundance than the members of the triad.

3.3.3 An example of an anti-synergistic effect

A sixth spot (367) shows strong induction by lovastatin (two-to threefold), and about half as much induction with lovastatin plus cholestyramine, but without sharing the animal-animal heterogeneity pattern of the 235-set (Fig. 13). This protein is also mitochondrial, and represents the clearest example of an anti-synergistic effect of lovastatin and cholestyramine. The existence of such an effect demonstrates that lovastatin and cholestyramine do not act exclusively through the same regulatory pathway.

3.3.4 Complexity of the cholesterol synthesis pathway

Taken together, these results suggest that treatment with lovastatin alone can affect both cytosolic and mitochondrial pathways using HMG-CoA, while cholestyramine, on the other hand, either alone or in combination with lovastatin, produces a strong effect on the putative cytosolic pathway, but little or no effect on the putative mitochondrial pathway. An explanation for this difference may lie in lovastatin's effect on levels of HMG-CoA and related precursor compounds that are exchanged between the cytosol and the mitochondrion, whereas cholestyramine should affect only the cytosolic pathways directly controlled by cholesterol and bile acid levels. It remains to be explained why some

proteins of the putative mitochondrial pathway are so much more variable in their expression in all groups. An examination of all the coregulated groups suggests that quantitative statistical techniques can extract a wealth of interesting information from large sets of reproducible gels. The abundance of spots in the 413 coregulation group, for example, shows an amazing level of concordance in their relative expression among the five individuals of the lovastatin and cholestyramine treatment group. This effect is not due to differences in total protein loading, since they have already been removed by scaling, and since proteins with quite different regulation patterns can be demonstrated (e.g., Fig. 13). Such effects raise the possibility that many gene coregulation sets may be revealed through the study of a sufficiently large population of control animals (i.e., without any experimental manipulation). This approach, exploiting natural biological variation in protein expression instead of drug effects, offers an important incentive for the construction of a large library of control animal patterns.

4 Conclusions

Because of the widespread use of rat liver in both basic biochemistry and in toxicology, there is a long-term need for a comprehensive database of liver proteins. The rat liver master pattern presented here has proven to be an accurate representation of this system, having been matched to more than 700 gels to date. As the number of proteins identified and the number of compounds tested for gene expression effects grows, we expect this database to contribute valuable insights into gene regulation. Its practical utility in several areas of mechanistic toxicology is already being demonstrated.

Received September 11, 1991

5 References

- [1] O'Farrell, P., J. Biol. Chem. 1975, 250, 4007-4021.
- [2] Klose, J., Humangenetik 1975, 26, 231-243.
- [3] Scheele, G. A., J. Biol. Chem. 1975, 250, 5375-5385.
- [4] Iborra, G. and Buhler, J. M., Anal. Biochem. 1976, 74, 503-511.
- [5] Anderson, N. G. and Anderson, N. L., Behring. Inst. Mitt. 1979, 63, 169-210.
- [6] Anderson, N. G. and Anderson, N. L., Clin. Chem. 1982, 28, 739-748.
- [7] Heydorn, W. E., Creed, G. J. and Jacobowitz, D. M., J. Pharmacol. Exp. Therap. 1984, 229, 622-628.
- [8] Anderson, N.L., Nance, S.L., Tollaksen, S.L., Giere, F.A. and Anderson, N. G., Electrophoresis 1985, 6, 592-599.
- [9] Racine, R. R. and Langley, C. H., Biochem. Genet. 1980, 18, 185-197.
- [10] Klose, J., Mol. Evol. 1982, 18, 315-328.
- [11] Neel, J. V., Baier, L., Hanash, S. and Erickson, R. P., J. Hered. 1985, 76, 314-320.
- [12] Marshall, R. R., Raj, A. S., Grant, F. J. and Heddle, J. A., Can. J. Genet. Cytol. 1983, 25, 457-446.
- Taylor, J., Anderson, N. L., Anderson, N. G., Gemmell, A., Giometti, C. S., Nance, S. L. and Tollaksen, S. L., in: Dunn, M. J. (Ed.), Electrophoresis '86, Verlag Chemie, Weinheim 1986, pp. 583-587.
- [14] Giometti, C. S., Gemmell, M. A., Nance, S. L., Tollaksen, S. L. and Taylor, J., J. Biol. Chem. 1987, 262, 12764-12767.
- [15] Anderson, N. L., Giere, F. A., Nance, S. L., Gemmell, M. A., Tollaksen, S. L. and Anderson, N. G., in: Galteau, M.-M. and Siest, G. (Eds.), Progrès Récents en Electrophorèse Bidimensionelle, Presses Universitaires de Nancy, Nancy 1986, pp. 253-260.
- [16] Anderson, N. L., Swanson, M., Giere, F. A., Tollaksen, S., Gemmell, A., Nance, S. L. and Anderson, N. G., Electrophoresis 1986, 7, 44-48.

- [17] Anderson, N. L., Giere, F. A., Nance, S. L., Gemmell, M. A., Tollaksen, S. L. and Anderson, N. G., Fundam. Appl. Toxicol. 1987, 8,39-50.
- [18] Anderson, N. L., in: New Horizons in Toxicology, Eli Lilly Symposium, 1991, in press.
- [19] Antoine, B., Rahimi-Pour, A., Siest, G., Magdalou, J. and Galteau, M. M., Cell. Biochem. Funct. 1987, 5, 217-231.
- [20] Elliott, B. M., Ramasamy, R., Stonard, M. D. and Spragg, S. P. Biochim. Biophys. Acta 1986, 870, 135-140.
- [21] Huber, B. E., Heilman, C. A., Wirth, P. J., Miller, M. J. and Thorgeirsson, S. S., Hepatology 1986, 6, 209-219.
- [22] Wirth, P. J. and Vesterberg, O., Electrophoresis 1988, 9, 47-53.
- [23] Witzmann, F. A. and Parker, D. N., Toxicol. Lett. 1991, 57, 29-36.
- [24] Rampersaud, A., Waxman, D. J., Ryan, D. E., Levin, W. and Walz, F. G., Jr., Arch. Biochem. Biophys. 1985, 243, 174-183.
- [25] Vlasuk, G. P. and Walz, F. G., Jr., Anal. Biochem. 1980, 105, 112-120.
- [26] Anderson, N. G. and Anderson, N. L., Anal. Biochem. 1978, 85, 331-
- [27] Anderson, N. L. and Anderson, N. G., Anal. Biochem. 1978, 85, 341— 354.
- [28] Anderson, L., Hofmann, J.-P., Anderson, E., Walker, B. and Anderson, N. G., in: Endler, A. T. and Hanash, S. (Eds.), Two-Dimensional Electrophoresis, VCH Verlagsgesellschaft, Weinheim 1989, pp. 288-207
- [29] Anderson, L., Two-Dimensional Electrophoresis: Operation of the ISO-DALT[®] System. Large Scale Biology Press, Washington, DC 1988, ISBN 0-945532-00-8, 170pp.
- [30] Neuhoff, V., Stamm, R. and Eibl, H., Electrophoresis 1985, 6, 427-448.

- [31] Neuhoff, V., Arold, N., Taube, D. and Ehrhardt, W., Electrophoresis 1988, 9, 255-262.
- [32] Anderson, N. L. and Hickman, B. J., Anal. Biochem. 1979, 93, 312–320.
- [33] Sidman, K. E., George, D. E., Barker, W. C. and Hunt, L. T., Nucl. Acids Res. 1988, 16, 1869-1871.
- [34] Taylor, J., Anderson, N. L., Coulter, B. P., Scandora, A. E. and Anderson, N. G., in: Radola, B. J. (Ed.), Electrophoresis '79, de Gruyter, Berlin 1980, pp. 329-339.
- [35] Taylor, J., Anderson, N. L. and Anderson, N. G., in: Allen, R. C. and Arnaud, P. (Eds.), Electrophoresis '81, de Gruyter, Berlin 1981, pp. 383-400.
- [36] Anderson, N. L., Taylor, J., Scancora, A. E., Coulter, B. P. and Anderson, N. G., Clin. Chem. 1981, 27, 1807-1820.
- [37] Taylor, J., Anderson, N. L., Scandora, A. E., Jr., Willard, K. E. and Anderson, N. G., Clin. Chem. 1982, 28, 861-866.
- [38] Taylor, J., Anderson, N. L. and Anderson, N. G., Electrophoresis 1983, 4, 338-345.
- [39] Anderson, N. L. and Taylor, J., in: Proceedings of the Fourth Annual Conference and Exposition of the National Computer Graphics Association, Chicago, June 26–30, 1983, pp. 69–76.
- [40] Anderson, N. L., Hofmann, J.-P., Gemmell, A. and Taylor, J., Clin. Chem. 1984, 30, 2031-2036.
- [41] Anderson, L., in: Schafer-Nielsen, C. (Ed.), Electrophoresis '88, VCH Verlagsgesellschaft, Weinheim 1988, pp. 313-321.
- [42] Neidhardt, F. C., Appleby, D. A., Sankar, P., Hutton, M. E. and Phillips, T. A., Electrophoresis 1989, 10, 116-121.
- [43] Gil, G., Goldstein, J. L., Slaughter, C. A. and Brown, M. S., J. Biol. Chem. 1986, 261, 3710-3716.

6 Addendum 1: Figures 1-13



Figure 1. Synthetic representation of the standard rat liver 2-D master pattern, rendered as a greyscale image using a videoprinter.

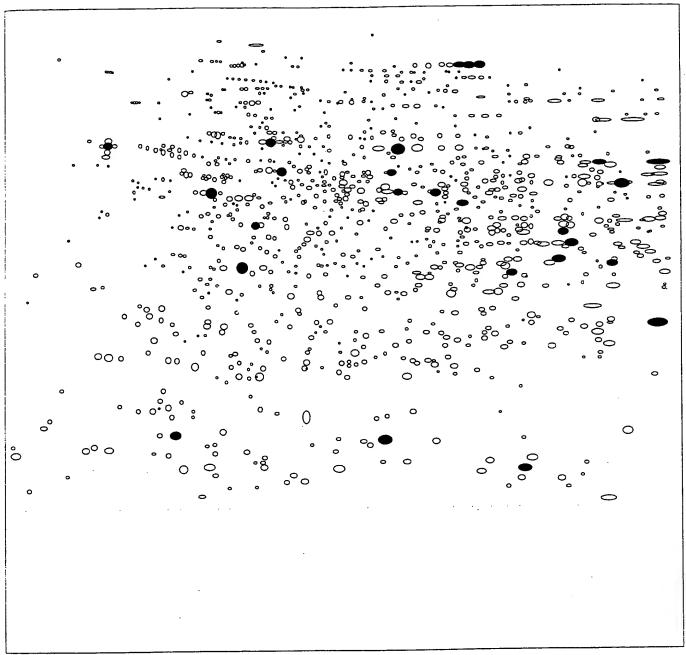


Figure 2. Schematic representation of the master pattern (the same as Fig. 1), useful as an aid in relating specific areas of Fig. 1 and the following detailed quadrants.

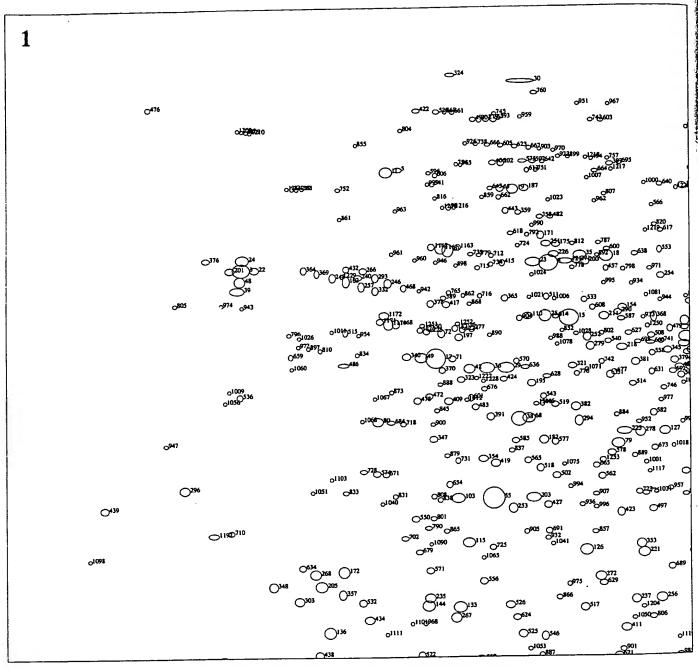


Figure 3. Upper left (high molecular weight, acidic) quadrant (#1) of the rat liver map, showing spot numbers.

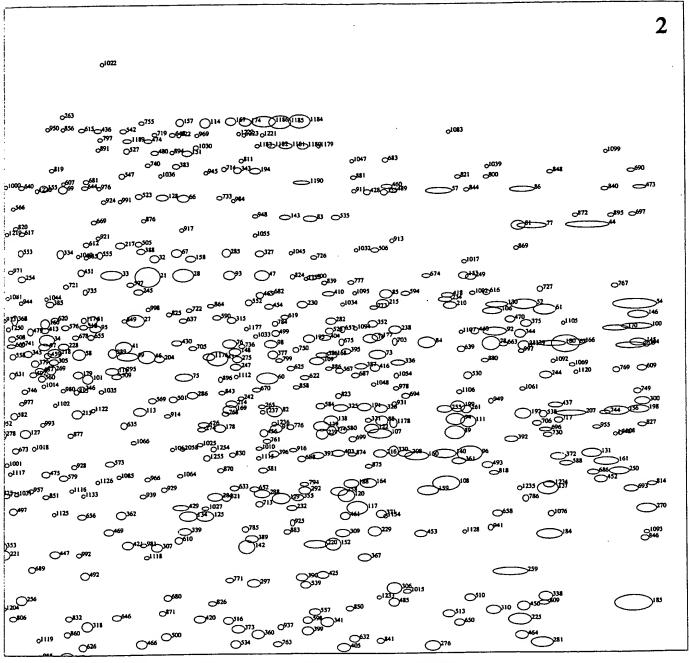


Figure 4. Upper right (high molecular weight, basic) quadrant (#2) of the rat liver map, showing spot numbers.

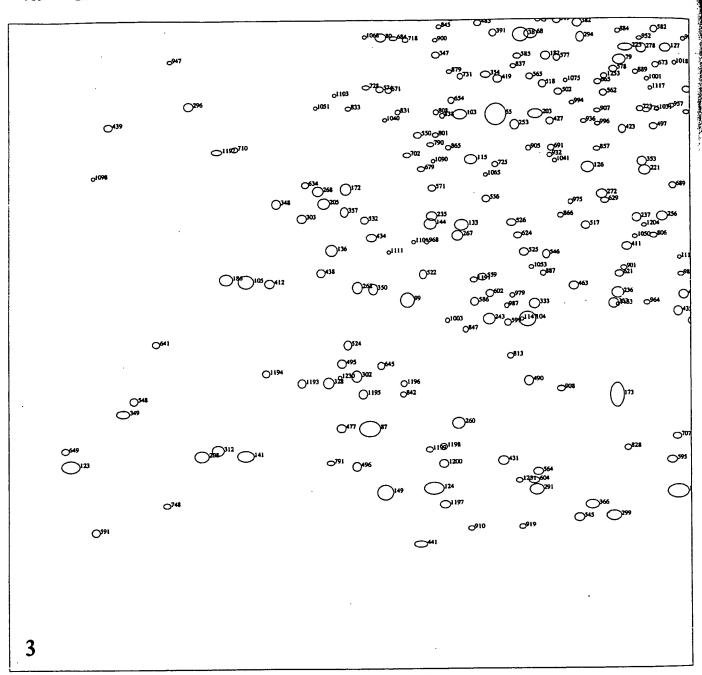


Figure 5. Lower left (low molecular weight, acidic) quadrant (#3) of the rat liver map, showing spot numbers.

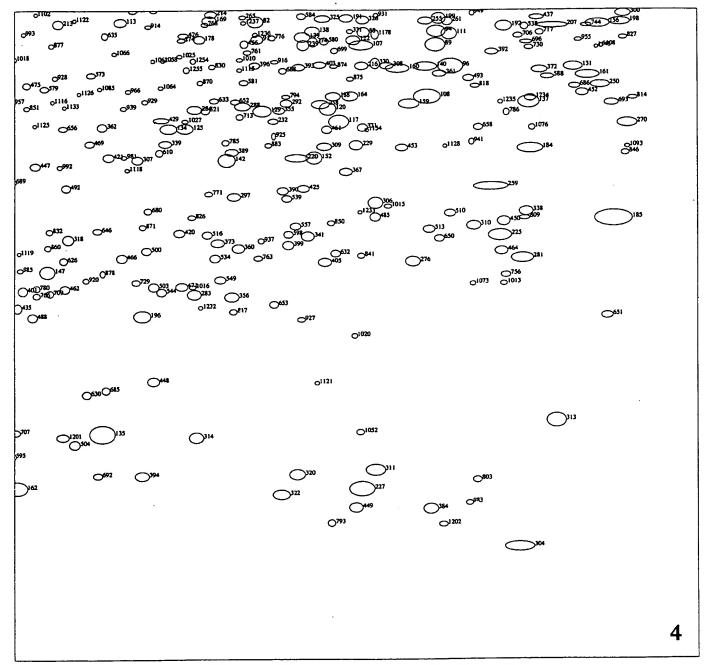


Figure 6. Lower right (low molecular weight, basic) quadrant (#4) of the rat liver map, showing spot numbers.

1250

750

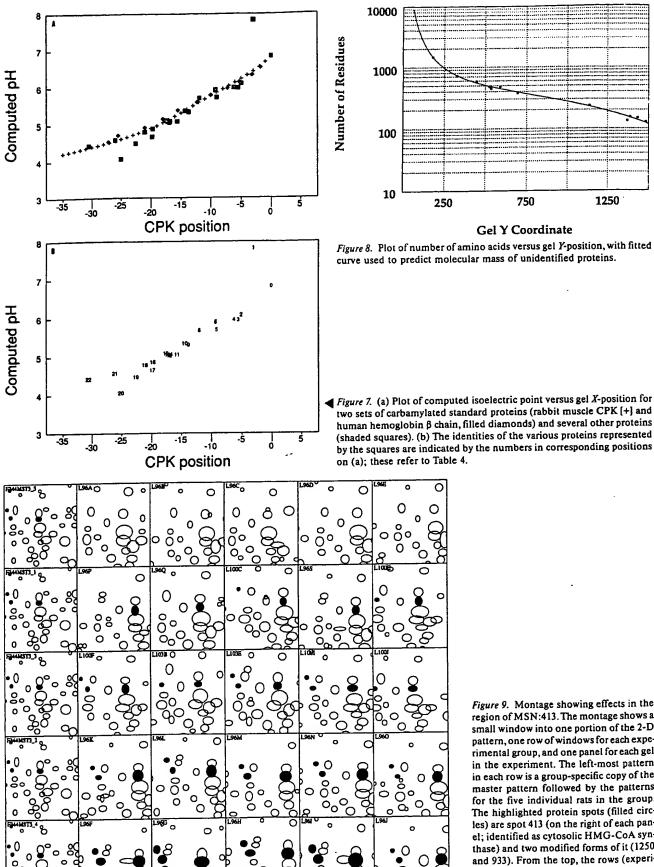


Figure 9. Montage showing effects in the region of MSN:413. The montage shows a small window into one portion of the 2-D pattern, one row of windows for each experimental group, and one panel for each gel in the experiment. The left-most pattern in each row is a group-specific copy of the master pattern followed by the patterns for the five individual rats in the group. The highlighted protein spots (filled circles) are spot 413 (on the right of each panel; identified as cytosolic HMG-CoA synthase) and two modified forms of it (1250 and 933). From the top, the rows (experimental groups) are: high cholesterol, controls, cholestyramine, lovastatin, and lovastatin plus cholestyramine.

Regulation of Rat Liver 413

(Putative Cytosolic HMG-CoA Synthase, 53kd) Test Compounds in Diet

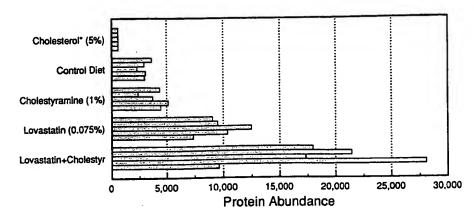


Figure 10. Bargraph showing the quantitative effects of various treatments on the abundance of MSN:413 (cytosolic HMG-CoA synthase) in the gels of Fig. 9.

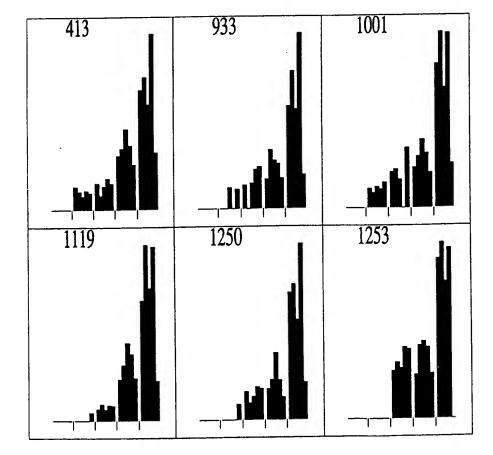


Figure 11. Bargraphs of a series of six coregulated spots including MSN:413. In the bargraphs, the abundances of the appropriate spot (master spot number shown at the top of the panel) in each animal are shown. The five five-animal groups are in the order (left to right): high cholesterol, controls, cholestyramine, lovastatin, and lovastatin plus cholestyramine. Each bar within a group represents one experimental animal liver (one 2-D gel). Note the correlated expression of the 6 spots, especially in the two far right (most strongly induced) groups.

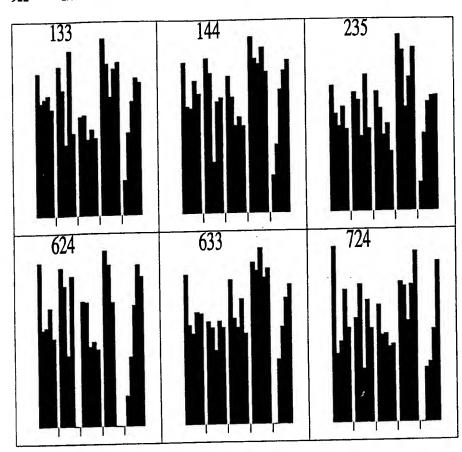


Figure 12. Data on a second coregulated group of spots, presented as in Fig. 11. The fourth experimental group (lovastatin) shows a modest induction, while the fifth group (lovastatin plus cholestyramine) does not.

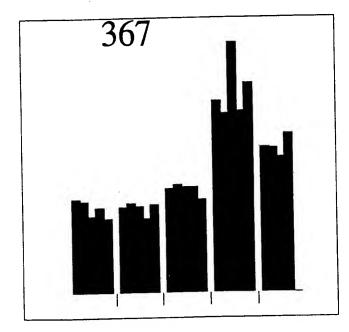


Figure 13. Data on spot MSN:367, presented as in Fig. 11. This protein shows unambiguously the anti-synergistic effect of lovastatin and cholestyramine (fifth group) as compared to lovastatin (fourth group). This response contrasts strongly with the regulation pattern seen in Fig. 11.

			proteins	000171	2.000	×	Y	СРКо	SDSMW	MSN	X	Y	CPKol	SDSMV
SN —	X	Y	CPKol	SDSMW	MSN	X	<u> </u>							
3	311	434	<-35.0	63,800	95	1119	536	-9.9	53,800	174	1364	183	-6.7	162,90
5	568	263	-24.3	102,900	96	1731	756 560	-2.0	40,700 51,600	175 177	825 1582	393 553	-15.7 -3.6	69,30 52,60
8	812	426	-16.0 <i>-</i> 25.2	64,800 101,000	97 98	1033 1406	566 565	-11.4 -6.1	51,700	178	1321	710	-7.2	43,00
1 5	549 845	268 520	-25.2 -15.3	55,200	99	578	1149	-23.8	25,000	179	1089	615	-10.4	48,30
7	629	589	-21.6	50,000	100	2004	538	>0.0	53,700	180	1866	567	-0.5	51,60
8	906	414	-14.0	66,300	101	1106	623	-10.1	47,900	181	411	295	-32.1	91,20
9	755	298	-17.5	90,200	102	482	455	-28.5	61,300	182	804	730	-16.2	42,00
20	649	403	-20.9	67,900	103	665	830	-20.2	37,300	184 185	1860 1997	896 1017	-0.6 >0.0	34,50 29,80
21	1204	448	-8.7	62,100	104 105	773 312	1182 1117	-17.0 <-35.0	23,800 26,100	186	279	1113	<-35.0	26,3
22	332 787	434 424	<-35.0 -16.6	63,800 65,000	105	1769	509	-1.5	56,100	187	773	296	-17.0	90,8
23 24	313	417	<-35.0	66,000	107	1585	720	-3.6	42,500	188	1538	807	-4.2	38,4
25	807	516	-16.1	55,500	108	1692	807	-2.4	38,300	191	1560	674	-3.9	44,9
27	1184	524	-9 .0	54,900	109	1482	593	-4.8	49,700	192	1818	687	-0.9	44,2
28	1263	446	-8.0	62,400	110	778	516	-16.9	55,500	193	1469 1380	555 266	-5.0 -6.4	52,44 101,6
29	743	605	-17.8	49,000	111 113	1728 1191	700 680	-2.0 -8.9	43,500 44,500	194 195	784	632	-16.7	47,3
30	768	112 417	-17.2 -8.6	348,600 66,000	114	1298	185	-7. 5	160,800	196	1227	1185	-8.4	23,7
32 33	1216 11 4 5	445	-9.5	62,500	115	682	907	-19.6	34,100	197	667	553	-20.1	52,6
~ 34	1037	555	-11.3	52,400	116	1146	610	-9.5	48,700	198	2006	681	>0.0	44,5
35	863	412	-14.9	66,600	117	1548	849	-4.1	36,500	199	1711	674	-2.2	44,9
36	712	606	-18.7	48,900	118	1050	577	-11.1	50,800	200	872	424	-14.7	65,0
38	763	694	-17.3	43,800	120	1530	828	-4.3	37,400	201	292	435 253	<-35.0 -18.0	63,7 107,8
39	304	470	<-35.0	59,800	121 122	838 1572	423 712	-15.4 -3.8	65,200 42,900	202 203	736 786	829	-16.7	37,4
11 12	1165 684	569 607	-9.2 -19.6	51,400 48,800	123	23	1433	<-35.0	15,300	204	1224	589	-8.5	50,0
13	1318	589	-7.3	50,000	124	621	1474	-21.9	13,900	205	439	983	-30.9	31,1
14	1924	362	-0.1	74,600	125	1298	862	-7.5	36,000	206	1994	571	>0.0	51,3
16	1203	586	-8.7	50,200	126	872	921	-14.7	33,50C	207	1895	687	-0.3	44,2
17	1391	447	-6.3	62,300	127	1000	717	-12.0	42,600	208	240	1418	<-35.0	15,8
18	309	454	<-35.0	61,500	128	1229	311	-8.4	86,100 37,300	210 211	1700 902	499 517	-2.3 -14.1	57,0 55,4
19	605	587 505	-22.5	50,100 53,000	129 130	1422 1776	832 499	-5.8 -1.4	57,000	213	1087	684	-10.4	44,4
50 51	621 1113	535 522	-21.8 -10.0	53,900 55,000	131	1930	757	-0.1	40,700	214	1340	668	-7.0	45,2
52	1820	499	-0.9	57,000	132	660	537	-20.4	53,800	215	1591	495	-3.5	57,3
53	725	177	-18.3	170,800	133	666	1019	-20.2	29,700	216	1585	755	-3.6	40,7
54	2001	500	>0.0	56,900	134	1271	862	-7.9	36,000	217	1159	393	-9.3	69,3
55	722	830	-18.4	37,300	135	1161	1389	-9.3 ~~ 7	16,800	218	931 713	572 177	-13.5 -18.7	51,2 170,5
56	678	533	-19.8	54,100 89,000	136 137	453 1858	1063 823	-29.7 -0.6	28,100 37,700	219 220	1479	911	-4.9	33,9
57 58	1682 1091	302 580	-2.5 -10.3	50,600	138	1504	697	-4.6	43,700	221	965	927	-12.8	33,3
59	1171	585	-9.2	50,300	139	1488	707	-4.8	43,200	223	934	716	-13.5	42,7
50	1400	624	-6.2	47,800	140	1689	756	-2.4	40,700	225	1812	1045	-1.0	28,8
51	1853	508	-0.6	56,200	141	311	1417	<-35.0	15,800	226	821	411	-15.8	66,8
52	1888	567	-0.4	51,500	142	1366	915	-6.7	33,800	227	1586	1483 567	-3.6 -10.8	13,6 51,6
55	735	297	-18.1	90,500	143	1429 615	346 1017	-5.7 -22.1	77,900 29,800	228 229	1065 1577	890	-3.7	34,8
56 57	1263	312 407	-8.0 -8.1	85,900 67,300	144 145	2006	566	>0.0	51,600	230	1458	496	-5.2	57,3
57 58	1252 779	692	-16.8	43,900	146	2006	518	>0.0	55,300	232	1440	849	-5.5	36,5
5 9	1064	296	-10.8	90,800	147	1070	1108	-10.7	26,500	234	1692	489	-2.4	57,9
71	656	589	-20.6	50,000	148	1347	578	-6.9	50,800	235	618	1004	-22.0	30,3
72	638	545	-21.2	53,100	149	541	1481	-25.7	13,700	236	920	1138	-13.7	25,4
73	1582	583	-3.6	50,400	150	1645 1269	760 236	-2.8 -7.9	40,500 117,000	237 238	952 1611	1008 541	-13.1 -3.2	30,2 53,5
14	1570	556	-3.8	52,300 48,000	151 152	1507	911	-7.9 -4.5	33,900	239	1489	720	-4.8	42,5
75 76	1264 1338	621 564	-8.0 -7.0	51,800	153	1722	448	-2.1	62,100	240	501	448	-27.7	62,1
77	1833	363	-0.8	74,400	154	932	503	-13.5	56,600	241	1820	569	-0.9	51,4
78	1767	565	-1.5	51,700	155	1031	294	-11.4	91,400	242	1357	658	-6.8	45,8
79	925	738	-13.6	41,600	156	1970	684	>0.0	44,400	243	711	1182	-18.7	23,8
30	534	698	-26.1	43,600	157	1258	183	-8.1	162,400	244	1855	621	-0.6	48,0
31	1811	363	-1.0	74,500	158	1275	417 820	-7.8 -2.6	65,900 37,800	245 246	1189 551	474 459	-8.9 `-25.1	59,3 61,0
32	1412	681	-6.0 -5.0	44,500 77,500	159 160	1663 1034	527	-2.6 -11.4	37,800 54,600	240	1348	604	-25.1 -6.9	49,1
33 34	1471 1662	347 563	-5.0 -2.7	77,500 51,800	161	1953	771	>0.0	40,000	248	460	448	-29.3	62,1
35 35	1596	479	-2.7 -3.4	58,900	162	1020	1482	-11.6	13,700	249	1733	451	-1.9	61,8
36	1817	301	-0.9	89,100	164	1566	806	-3.8	38,400	250	1974	788	>0.0	39,2
B7	516	1371	-27.0	17,400	166	1905	565	-0.2	51,700	251	808	392	-16.1	69,5
88	1589	698	-3.5	43,600	167	1340	181	-7.0	164,900	252	874	553	-14.6	52,5
89	1706	719	-2.2	42,500	168	1506	583	-4.6	50,400	253	753	848	-17.6	36,5
90	651	329	-20.8	81,700	169	1338	678	-7.0 - 0.0	44,700	254 255	995 1690	450 679	-12.1 -2.4	61,9 44,6
91	1415	710	-6.0	43,000 53,200	170	1969	541 378	>0.0 -16.3	53,500 71,800	255 256	994	1006	-2.4 -12.1	30,2
92	1773 1338	545 446	-1.4 -7.0	53,200 62,300	171 172	800 476	958	-16.3 <i>-</i> 28.7	71,800 32,100	257	508	464	-12.1	50,2 60,4
93						710		-a-U./	Ja., 100			,		, -

^{a)} Master table of proteins in the rat liver database, showing spot master number, gel position (x and y), isoelectric point relative to CPK standards, and predicted molecular mass (from the standard curve of Fig. 8).

									<u> </u>					
MSN	x	Y	CPKol	SDSMW	MSN	X	Υ	СРКФІ	SDSMW	MSN	x	Y	CPKpl	SDSMW
250	1796	961	-1.1	31,900	345	1006	578	-11.9	50,800	426	1296	704	-7. 6	43,300
259 260	661	1361	-20.4	17,700	346	1095	640	-10.3	46,800	427 428	810 1565	843 303	-16.0 -3.9	36,800 88,700
261	1725	679	-2.0	44,600	347	625 361	728 983	-21.7 -35.3	42,000 31,100	429	1259	847	-8.0	36,600
262	496	1127	-28.0 -10.9	25,800 177,400	348 349	110	1343	<-35.0	18,300	430	1253	562	-8.1	51,900
263 265	1063 1390	172 673	-10.9	45,000	350	521	1130	-26.7	25,700	431	734	1426	-18.1 29.5	15,500
266	510	437	-27.3	63,400	351	912	619	-13.9	48,100 54,300	432 434	483 518	433 1041	-28.5 -26.9	63,900 28,900
267	660	1038	-20.4	29,000	352 353	1574 961	530 912	-3.7 -12. 9	33,900	435	1020	1170	-11.6	24,300
268	430	961 606	-31.0 -11.2	31,900 48,900	354	706	762	-18.9	40,400	436	1122	196	-9.8	147,600
269 270	1044 2019	853	>0.0	36,300	355	1450	830	-5.3	37,300	437	1870	673	-0.5	45,000 26,700
271	857	422	-15.0	65,200	356	1374	1152	-6.5	24,900 30,600	438 439	435 86	1102 847	-31.0 <-35.0	36,600
272	895	968	-14.2	31,700	357 358	474 798	997 346	-28.7 -16.3	77,800	440	1740	544	-1.8	53,200
274	1292	712 590	-7.6 -6.9	42,900 49,900	359	764	338	-17.3	79,400	441	599	1571	-22.8	10,800
275 276	1350 1670	1089	-2.6	27,100	360	1384	1068	-6.4	27,900	443	743	335	-17.8	80,100
277	688	538	-19.4	53,700	361	1713	769	-2.1	40,100	446 447	801 1050	668 926	-16.2 -11.1	45,200 33,300
278	961	718	-13.0	42,600	362	1161 914	859 1156	-9.3 -13.8	36,100 24,800	448	1245	1298	-8.2	19,800
279	879	570	-14.5 -0.7	51,300 27,300	363 364	412	435	-32.0	63,700	449	1576	1516	-3.7	12,600
281 282	1848 1505	1084 525	-4.6	54,800	365	741	486	-17.9	58,200	450	1818	1021	-0.9	29,600
283	1313	1147	-7.3	25,100	366	878	1503	-14.6	13,000	451	1094	440	-10.3 >0.0	63,100 38,600
284	1314	829	-7.3	37,400	367	1560	935	-3.9	33,000 55,200	452 453	1945 1652	802 894	-2.8	34,600
285	1332	408	-7.1	67,200	368 369	983 434	520 441	-12.4 -31.0	63,000	454	1403	500	-6.1	56,900
286	1277 1391	652 824	-7.8 -6.3	46,100 37,600	370	639	610	-21.2	48,700	456	1394	718	-6.3	42,600
288 289	1147	579	-9.5	50,700	371	1587	860	-3.6	36,100	457	905	436	-14.0	63,500
290	925	511	-13.6	55,900	372	1875	762	-0.5	40,400	459 460	1038 1598	581 294	-11.3 -3.4	50,500 91,400
291	787	1476	-16.6	13,900	373	1351	1059 715	-6.8 -4.6	28,300 42,700	461	1528	863	-4.3	35,900
292	1462	818	-5.1 -26.3	37,800 62,000	374 375	1506 1823	532	-0.9	54,200	462	1098	1137	-10.2	25,400
293 294	531 860	449 698	-20.3 -14.9	43,600	376	254	417	<-35.0	65,900	463	849	1125	-15.2	25,800
295	1162	609	-9 .3	48,700	377	1409	583	-6.1	50,400	464	1814	1072	-0.9	27,800 58,700
296	218	814	<-35.0	38,000	378	621	494	-21.8	57,500 49,600	465 466	1388 1194	481 1084	-6.3 -8.9	27,300
297	1377	979	-6.5 40.0	31,300 12,400	379 381	1017 953	595 598	-11.7 -13.1	49,400	468	577	467	-23.9	60,100
299 300	913 2012	1523 667	-13.9 >0.0	45,300	382	856	674	-15.0	44,900	469	1140	888	-9.6	34,900
301	702	178	-19.0	169,200	383	1252	258	-8.1	105,300	470	1797	524	-1.1 7.6	54,800 25,500
302	494	1280	-28.1	20,400	384	1699	1518	-2.3	12,500 57,500	471 472	1293 618	1133 655	-7.6 -21.9	46,000
303	403	1008	-32.6	30,100 10,300	385 386	1042 1490	493 583	-11.2 -4.7	50,400	473	2009	299	>0.0	89,900
304 305	1843 1049	1585 593	-0.7 -11.1	49,800	387	1554	603	-4.0	49,100	474	1205	215	-8.7	131,300
306	1608	989	-3.3	30,900	388	1193	404	-8.9	67,700	475	1035	788	-11.4	39,200 207,600
307	1219	916	-8.5	33,700	389	1374	902	-6.5	34,300	476 477	160 469	155 1370	<-35.0 -28.9	17,400
308	1627	755	-3.0	40,700	390 391	1456 718	969 690	-5.2 -18.5	31,700 44,000	478	599	662	-22.8	45,600
309 310	1524 1769	892 1028	-4.4 -1.5	34,700 29,400	392	1799	732	-1.1	41,900	479	1009	540	-11.8	53,500
311	1609	1451	-3.3	14,700	393	1482	758	-4.8	40,600	480	1216	235	-8.6	117,400
312	266	1408	<-35.0	16,100	394	1227	1461	-8.4	14,400	482	816	346 673	-15.9 -19.3	77,800 44,900
313	1902	1365	-0.3	17,600	395	1530 1410	577 755	-4.3 -6.0	50,800 40,800	483 485	693 1608	673 1013	-3.3	30,000
314	1316	1395 523	-7.3 -7.0	16,600 54,900	396 397	912	256	-13.9	106,400	486	478	599	-28.6	49,300
315 318	1341 1104	1053	-10.1	28,500	399	1465	1063	-5.0	28,100	487	1025	607	-11.5	48,800
320	1480	1459		14,400	400	1473	450	-4.9	61,900	488	1045	1186	-11.2	23,700 89,200
321	850	603		49,100	401	1029	1140	-11.5	25,300 40,800	489 490	1609 775	301 1289	-3.3 -17.0	20,100
322	1454	1494		13,300 47,700	403 404	1516 1495	754 554	-4.4 -4.7	52,500	491	692	178	-19.3	169,300
323	670 655	626 101		420,500	405	1525	1092		27,100	492	1100	964	-10.2	31,800
324 325	1521	675		44,800	406	723	252		108,000	493	1760	776		39,700
326	1587	677		44,700	409	650	663		45,500	494	882	247 1258		110,700 21,200
327	1388	409		67,000	410	1501	478		59,000 28,300	495 496	470 494	1436		15,200
328	448	1291	-30.0	20,100 40,900	411 412	936 350	1057 1120	-13.4 -35.9	26,000 26,000	497	980	852		36,400
330 331	1608 1566	751 697	-3.3 -3.8	43,700	413		538		53,700	499	1414	546	-6.0	53,100
332	531	471		59,600	415	737	425	-18.0	64,900	500	1234	1072		27,800 45,700
333	784	1156	-16.7	24,700	416		606		48,900 57,300	501 502	1246 824	659 792		45,700 39,000
334	1059	407		67,300	417		496 482		57,300 58,600	502 503	1246	1134		25,500
335	1593	303 598		88,500 49,400	418 419				40,000	504	1115	1407	-9.9	16,200
336 338	1616 1854	1004		30,300	420		1041		28,900	505	1189	391		69,700
339	1265	888		34,900	421	1171	912	-9.1	33,900	506	1578	402		68,000
340	581	585	-23.6		422				193,700	507 508	787 979	250 552		109,000 52,600
341	1497			28,700	423				36,200 47,700	508 509	1153	619		48,100
343	1351	265			424 425				31.800	510		1006		30,200
344	1813	549	-0.9	32,000	720	50		••••		_				

			0014-1	CDCMW	MSN	×	Υ	CPKpl	SDSMW	MSN	×	Y	CPKpl	SDSMW
MSN	×	Y	CPKoi	SDSMW	MON									
511	809	484	-16.0	58,400	596	619	269	-21.9	100,500	674	1661	448	-2.7	62,100
512	1099	533	-10.2	54,100	597	1176	461	-9.1	60,700	675	1523	562	-4.4 10.0	51,900 46,700
513	1696	1034	-2.3	29,200	598	1465	1044	-5.0	28,800	676 677	708 919	642 615	-18.8 -13.7	48,300
514	948	636	-13.2	47,100	599	741	1188	-17.9	23,600	677 678	1085	551	-13.7 -10.5	52,700
515	481	543	-28.5	53,400	600	907	402	-14.0	68,000 45 BOO	679	600	923	-22.7	33,400
516	1334	1044	-7.1	28,800	601	687	658	-19.5	45,800 25,400	680	1237	1004	-8.3	30,300
517	868	1021	-14.8	29,700	602	712	1138	-18.7 -14.1	165,200	681	1103	283	-10.1	95,100
518	798	779	-16.3	39,600	603 604	898 783	181 1461	-16.7	14,400	682	1406	477	-6.1	59,100
519	822	670	-15.7 -21.5	45,100 189,000	605	736	223	-18.0	125,300	683	1596	249	-3.4	109,800
520	632	165 830	-21.5 -7.1	37,300	606	629	273	-21.6	98,700	684	555	699	-24.8	43,500
521 522	1332 603	1104	-22.6	26,600	607	1064	286	-10.8	94,000	685	1167	1313	-9.2	19,300
523	1190	309	-8.9	86,800	608	883	503	-14.5	56,700	686	1932	790	0.0	39,100
524	479	1226	-28.6	22,300	609	2012	610	>0.0	48,700	687	1545	619	-4.1	48,100
525	768	1066	-17.2	28,000	610	1255	903	-8.1	34,200	688	1456	764	-5.2	40,300 32,300
526	747	1016	-17.7	29,800	612	1103	391	-10.1	69,600	689	1011	953	-11.8	100,200
527	1170	231	-9 .2	119,600	613	778	265	-16.9	102,000	690	1995	270 888	>0.0 -16.0	34,900
528	1502	542	-4.6	53,400	614	*824	518	-15.7	55,400	691	812 1154	1461	-10.0	14,400
530	1728	620	-2.0	48,000	615	1095	195	-10.3	149,100	692 693	1993	819	>0.0	37,800
532	507	1011	-27.4	30,000	616	1759 994	478 372	-1.6 -12.1	59,000 72,900	694	1628	656	-3.0	45,900
533	870	489	14.7	57,900 07,300	617 618	751	374	-17.6	72,400	695	928	254	-13.6	107,000
534	1347	1085	-6.9	27,300	619	1429	518	-5.7	55,300	696	1854	715	-0.6	42,700
535	1513	346	-4.5 <-35.0	77,800 46,000	620	1050	520	-11.1	55,200	697	1997	345	>0.0	78,000
536 536	308 1851	654 689	-0.7	44,100	621	923	1105	-13.7	26,600	698	957	563	-13.0	51,800
538 539	1463	982	-5.1	31,100	622	1462	622	-5.1	47,900	699	1540	730	-4.2	42,000
540	909	561	-13.9	52,000	623	759	225	-17.4	124,000	702	577	900	-23.8	34,400
541	625	289	-21.7	93,100	624	758	1038	-17.4	29,000	703	1610	562	-3.2	51,900
542	1164	198	-9.2	146,200	625	1438	606	-5.5	48,900	705	1278	571	-7.8 -0.7	51,200 43,300
543	803	655	-16.2	45,900	626	1096	1089	-10.2	27,200	706 707	1841 1018	704 1386	-0.7 -11.7	16,900
544	1259	1143	-8.0	25,200	627	942	548	-13.3	53,000 48,000	707	1074	1145	-10.7	25,100
545	856	1526	-15.0	12,200	628	809	621 979	-16.0 -14.1	31,300	710	293	889	<-35.0	34,800
546	803	1071	-16.2	27,800	629 630	899 1135	1321	-9.6	19,100	712	720	412	-18.5	66,600
547	1162	274	-9.3 - 35.0	98,400 19,000	631	979	615	-12.5	48,300	713	1386	841	-6.4	36,800
548	128 1355	1321 1122	<-35.0 -6.8	25,900	632	1542	1076	-4.1	27,600	714	1328	263	-7.1	103,100
549 550	595	866	-23.0	35,800	633	1345	814	-6.9	38,000	715	698	433	-19.1	63,900
552	1369	494	-6.6	57,500	634	409	950	-32.2	32,400	716	701	481	-19.0	58,700
553	992	405	-12.2	67,600	635	1165	704	-9.2	43,300	717	1875	699	-0.5	43,600
555	1125	410	-9.8	66,900	636	774	604	-17.0	49,000	718	575	702	-23.9 -8.6	43,400 140,400
556	705	975	-18.9	31,400	637	1263	524	-8.0	54,800	719 721	1216 1069	204 464	-10.8	60,400
557	1477	1030	-4.9	29,300	638	952	411	-13.1	66,700 51,000	721	1272	506	-7.9	56,400
558	980	583	-12.5	50,400	639	1717	575 292	-2.1 -12.1	92,000	723	958	822	-13.0	37,700
559	700	1109	-19.1	26,400	640 641	· 994 165	1224	<-35.0	22,400	724	763	395	-17.3	69,100
560	1028	621	-11.5	48,000 38,900	642	803	251	-16.2	108,900	725	720	916	-18.5	33,700
562	898	794	-14.1 -16.6	14,900	643	719	296	-18.5	90,700	726	1476	415	-4.9	66,200
564	789 777	1446 766	-16.9	40,200	644	1100	294	-10.2	91,400	727	1846	473	-0.7	59,400
565 566	980	328	-12.5	81,900	645	534	1263	-26.1	21,000	728	510	783	-27.3	39,400
567	1519	611	-4.4	48,600	646	1153	1038	- 9 .4	29,000	729	1217	1126	-8.6	25,800
569	1212	661	-8.6	45,600	648	1246	204	-8.2	140,000	730	1858	724	-0.6	42,300 40,300
570	760	594	-17.4	49,700	649	14	1406		16,200	731	665	765 312	-20.2 -7.2	85,900
571	618	956	-21.9	32,100	650	1713	1049		28,600	733	1321 719	427	-18.5	64,600
573	1142	771	-9.6	40,000	651	1986	1183	>0.0	23,800	734 735	1101	473	-10.2	59,500
574	532	787	-26.2	39,300	652	1378	816		38,000 24,400	736	1359	569	-6.7	51,400
575	771	250	-17.1	109,200	653	1442	1165 806		38,400	738	696	220	-19.2	127,600
576	1068	534	-10.8	54,100	654 655	650 1111	551		52,700	739	687	409	-19.5	67,000
577	822	734	-15.7	41,800 40,800	656	1095	861		36,000	740	1205	256	-8.7	106,200
578	914	754 704	-13.8	38,900	657	1524	540		53,600	741	995	563	-12.1	51,900
579	1064	794	-10.8 -4.4	42,800	658	1777	860		36,000	742	898	596	-14.1	49,500
580	1524 1392	714 783		39,400	659	391	584		50,400	743	881	181	-14.5	165,900
581 582	982	686		44,200	660	977	565		51,700	744	1951	686	_	44,200
584	1487	672		45,000	661	658	166	-20.5	187,500	745	726	168		183,600
585	758	731	-17.4	41,900	662		312		86,100	746	999	643		46,600 13,000
586	687	1152	_	24,900	663		567		51,500	748	182	1503 649		46,300
587	930	523		55,000	664	888	268		100,900	749 760	2005 1448	575		51,000
588	1888	774	-0.4	39,900	665		775		39,800	750 751	792	266		101,900
589	642	485		58,300	666				126,300	751 752		296		90,600
590	1317	519		55,300	667		227		122,400 189,100	752 754	664	254		107,000
591	65	1548	_	11,500	668				76,300	755	1195	184		161,000
592	1014	614		48,400 172,300	669 670				46,600	756	1821	1113		26,300
593 504	732			172,300 59,000	671	547			39,200	757		246	-13.9	111,000
594 595	1627 1009	478 1426		15.500	673				41.200	760		133	-16.5	264,900
292	1003	1420	-11.0	, 5.566	5.0	301								

MSN	Х	Y	СРКЫ	SDSMW	MSN	X	Y	CPKpI	SDSMW	MSN	X	Y	CPKol	SDSMW
				44.000	040	1863	271	-0.6	99,500	939	1197	827	-8.8	37,500
761	1399	733	-6.2	41,800	848		523	-9.2	54,900	941	1765	885	-1.5	35,000
763	1416	1085	-5.9	27,300	849	1166		-4.2	29,600	942	602	472	-22.7	59,600
764	2020	569	>0.0	51,400	850	1535	1024	-11.4	37,500	943	312	498	<-35.0	57,100
765	651	475	-20.8	59,300	851	1035	826		53,400	944	993	491	-12.1	57,700
766	1052	1149	-11.1	25,000	852	834	542	-15.5		945	1300	269	-7.5	100,300
767	1968	468	>0.0	59,900	855	499	220	-27.8	127,100	946	630	423	-21.6	65,100
768	1330	685	-7.1	44,300	856	1063	194	-10.9	150,500				<-35.0	41,600
769	1970	613	>0.0	48,500	857	887	890	-14.4	34,800	947	187	736	-6.5	78,200
770	857	617	-15.0	48,200	858	1448	639	-5.4	46,900	948	1380	344		
	1337	974	-7.0	31,500	859	706	311	-18.9	86,200	949	1766	665	-1.5	45,400
771		502	-3.7	56,700	860	1070	1066	-10.7	28,000	950	1038	193	-11.3	151,000
773	1576	824	-12.8	37,600	861	472	347	-28.8	77,600	951	860	152	-14.9	213,000
775	969	708	-5.5	43,100	862	674	480	-19.9	58,800	952	957	701	-13.0	43,400
776	1438	458	-4.2	61,000	864	1307	499	-7.4	57,000	954	503	547	-27.6	53,000
777	1539	434	-15.1	63,800	865	645	887	-21.0	34,900	955	1938	712	>0.0	42,900
778	850	411	-19.1	66,800	866	827	1004	-15.6	30,300	957	1010	816	-11.8	37,900
779	700		-11.1	25,500	868	685	494	-19.5	57,400	959	768	174	-17.2	174,900
780	1052	1136		54,400	869	1807	402	-1.0	68,000	960	596	419	-23.0	65,700
784	1413	529	-6.0		870	1323	783	-7.2	39,400	961	557	409	-24.8	67,100
785	1364	885	-6.7	35,000	871	1228	1031	-8.4	29,300	962	887	320	-14.4	83,900
786	1822	835	-0.9	37,100	872	1904	346	-0.3	77,700	963	564	334	-24.5	80,500
787	893	392	-14.3	69,500		556	647	-24.8	46,400	964	969	1155	-12.8	24,800
790	616	882	-22.0	35,100	873			-4.2	40,700	965	671	255	-20.0	106,600
791	451	1429	-29.8	15,400	874	1540	756		39,700	966	1204	798	-8.7	38,700
792	777	377	-16.9	72,000	875	1566	777	-3.8	•	967	910	154	-13.9	210,300
793	1536	1543	-4.2	11,700	876	1198	351	-8.8	76,800	968	609	1048	-22.3	28,700
794	1461	807	-5.1	38,300	877	1076	720	-10.6	42,500		1285	206	-7.7	138,900
796	388	546	-33.6	53,100	878	1161	1111	-9.3	26,400	969 970	822	232	-15.8	119,300
797	1126	212	-9.8	133,700	879	647	757	-20.9	40,700			437	-12.6	63,400
798	933	437	-13.5	63,400	880	1756	594	-1.6	49,700	971	976		-32.6	51,600
799	1420	593	-5.9	49,800	881	1543	278	-4.1	97,100	972	403	567		57,400
800	1759	279	-1.6	96,500	883	1432	890	-5.7	34,800	974	279	495	<-35.0	31,200
801	624	865	-21.7	35,800	884	922	689	-13.7	44,100	975	844	981	-15.3	
802	898	547	-14.2	53,000	885	1103	414	-10.1	66,400	976	1124	295	-9.8	91,100
803	1775	1468	-1.4	14,200	886	1501	607	-4.6	48,900	977	994	664	-12.1	45,400
804	573	196	-24.0	148,400	887	798	1103	-16.3	26,600	978	1612	642	-3.2	46,700
805	203	494	<-35.0	57,400	888	636	634	-21.3	47,200	979	749	1141	-17.7	25,300
806	980	1039	-12.5	29,000	889	951	759	-13.1	40,600	980	1064	642	-10.8	46,700
807	902	308	-14.1	87,200	890	717	548	-18.6	52,900	981	1197	911	-8.8	33,900
808	625	827	-21.7	37,500	891	1123	229	-9.8	121,200	983	1762	1508	-1.6	12,800
		1015	-0.7	29,900	892	891	413	-14.3	66,400	984	1344	317	-6.9	84,700
809	1851	573	-30.9	51,100	894	1245	234	-8.2	117,800	985	1024	1105		26,600
810	440	249	-6.8	109,700	895	1962	346	>0.0	77,700	987	739	1159	-17.9	24,600
811	1358		-15.1	69,400	896	1322	626	-7.2	47,700	988	816	555	-15.9	52,400
812	851	393		21,600	897	420	570	-31.4	51,300	990	785	361	-16.7	74,900
813	745	1246	-17.8	38,200	898	662	428	-20.3	64,500	991	1159	317	-9.3	84,500
814	2028	810			899	845	243	-15.3	113,000	992	1090	928	-10.4	33,300
815	1086	645		46,500 05 700	900	624	703	-21.7	43,400	993	1030	701	-11.5	43,400
816	629	313	_	85,700	901	931	1094	-13.5	27,000	994	847	811	-15.2	38,200
817	1376	1177		24,000		799	229	-16.3	121,000	995	902	461	14.1	60,700
818	1771	790		39,100	903	765	520		55.200	996	888	847		36,600
819	1045	263		103,100	904				34,800	997	1815	579		50,700
820	984	362		74,600	905	775	889 824		37,600	998	1205	504		56,500
821	1712	279		96,700	907	888			19,700	999	617	289		93,100
822	1256	205		139,200	908	828	1303		11,700	1000	968	290		92,700
823	1517	654		46,000	910		1544		89,100	1001	970	771		40,000
824	1442	449		62,000	911	1544	301	-4.1	70,400	1002	1736	478		58,900
825		513		55,800	913	1606	387			1002	643	1184		23,700
826	1309	1014		29,900	914		688		44,100	1005	822	487		58,100
827	2012	708		43,100	916		749		41,100		875	279		96,400
828	937	1405	-13.4	16,200	917		367		73,700	1007	291	644		46,600
830	1342	756	-7.0	40,700	919		1541	_	11,700	1009				41,200
831	562	826	-24.5	37,500	920		1123		25,900	1010		745		53,500
832		1039		29,000	921		380		71,500	1011	459	541		45,600
833		820		37,800	923		242	_	113,200	1012		661		
834		581		50,500	924	1131	318		84,300	1013		1128		25,800 47,200
837		748			925	1441	874	-5.5	35,400	1014		634		47,200
838					926		. 219	-19.7	128,200	1015		994		30,700
839					927		1191		23,500	1016		1134		25,500
840					928		775		39,800	1017		424		
					929		816		38,000	1018		743		
841					931		670		45,100	1020	1574	1219		
842		1312			932				34,400	1021	781	484		
843					933				55,100	1022	1129	83		
844					934				60,600	1023		317		
845					936				36,800	1024		446	-16.7	62,400
846					937				28,400	1025		739	-7.7	41,500
847	673	1200	-19.9	23,200	937	1441	1030	3.5	_0,			-		

MSN	×	Υ	CPKol	SDSMW	MSN	×	Υ	CPKpI	SDSMW
									
1026	405	552	-32.5	52,600	1153	921	1158	-13.7	24,700
1027	1298	848	-7.5	36,500	1154	1594	864	-3.5	35,900
1028	856	547	-15.0	53,000	1161	637 623	400 397	-21.3	68,400 68,800
1030	1284	226	-7.7 12.2	123,200	1162 1163	665	397	-21.8 -20.2	68,700
1031	986	822	-12.3 -4.1	37,700 67,900	1168	564	528	-24.4	54,500
1032	1547 1381	403 551	-6.4	52,700	1170	552	529	-25.0	54,500
1033 1034	1525	496	-4.3	57,200	1171	538	524	-25.9	54,800
1035	1128	645	-9.7	46,500	1172	545	514	-25.5	55,700
1036	1226	274	-8.5	98,300	1174	1099	522	-10.2	55,000
1039	1761	262	-1.6	103,600	1176	1304	586	-7.5	50,200
1040	541	839	-25.7	36,900	1177	1366	539	-6.6	53,700
1041	818	910	-15.8	34,000	1178	1608	702	-3.3	43,400
1044	1036	485	-11.3	58,300	1179	1485	224	-4.8	124,900
1045	1439	407	-5.5	67,300	1180	1459	224	-5.2 5.7	124,900
1047	1540	250	4.2	109,200	1181	1431 1407	223 223	-5.7 -6.1	125,100 125,200
1048	1576	635	-3.7	47,100 se 700	1182 1183	1383	224	-6.4	124,700
1049	1089	411 1040	-10.4 -13.2	66,700 28,900	1184	1454	182	-5.3	164,400
1050 1051	949 426	818	-31.1	37,800	1185	1422	183	-5.8	162,600
1052	1583	1385	-3.6	16,900	1186	1394	182	-6.3	164,300
1053	779	1092	-16.8	27,000	1189	1171	214	-9.2	131,800
1054	1613	620	-3.2	48,000	1190	1457	286	-5.2	94,200
1055	1380	377	-6.5	72,000	1191	686	1114	-19.5	26,200
1056	284	663	<-35.0	45,500	1192	265	893	<-35.0	34,700
1058	1261	746	-8.0	41,200	1193	403	1292	-32.6	20,000
1060	393	605	-33.3	49,000	1194	344	1275	<-35.0	20,600
1061	1817	645	-0.9	46,600	1195	505 572	1311 1293	-27.6 -24.1	19,400 20,000
1062	1245	746	-8.2	41,200	1196 1197	639	1502	-21.2	13,000
1064	1258	792 934	-8.1 -18.9	39,000 33,000	1198	637	1402	-21.3	16,300
1065 1066	705 1181	734	-9.0	41,800	1199	614	1407	-22.1	16,200
1067	529	658	-26.3	45,800	1200	637	1431	-21.3	15,400
1068	508	696	-27.4	43,700	1201	1095	1394	-10.3	16,600
1069	1898	604	-0.3	49,100	1202	1719	1545	-2.1	11,600
1071	873	609	-14.7	48,700	1203	791	668	-16.5	45,200
1073	1768	1128	-1.5	25,800	1204	964	1021	-12.9	29,700
1075	836	773	-15.4	39,900	1205	313	195	<-35.0	148,700 149,800
1076	1863	861 566	-0.6	36,000 51,600	1208 1209	306 320	194 197	<-35.0 <-35.0	147,400
1078	826 971	566 483	-15.7 -12.7	51,600 58,500	1210	326	197	<-35.0	146,600
1081 1083	1697	202	-2.3	142,300	1211	394	294	-33.2	91,400
1085	1157	794	-9.4	38,900	1212	402	294	-32.7	91,200
1090	620	910	-21.9	34,000	1214	386	294	-33.7	91,400
1092	1867	597	-0.5	49,500	1215	641	329	-21.2	81,600
1093	2019	894	>0.0	34,600	1216	660	329	-20.4	81,600
1094	1546	538	-4.1	53,700	1217	914	266	-13.8	101,800
1095	1545	477	-4.1	59,100	1218	873	245	-14.7	112,000
1098	61	935	<-35.0	33,000	1219	970	372	-12.7	72,900
1099	1954	237	>0.0	116,000	1220	1021	298 205	-11.6 -6.3	90,100 139,500
1101	588	1048	-23.3	28,600	1221 1222	1392 1354	203	-6.8	141,800
1102 1103	1050 457	667 797	-11.1 -29.5	45,200 38,800	1223	1362	205	-6.7	139,500
1105	1884	532	-0.4	54,200	1224	673	540	-19.9	53,600
1106	1714	649	-2.1	46,300	1225	614	542	-22.1	53,400
1107	1717	546	-2.1	53,100	1226	603	539	-22.6	53,600
1108	1976	722	>0.0	42,400	1227	696	623	-19.2	47,800
1111	547	1066	-25.3	28,000	1228	707	628	-18.9	47,500
1112	1348	621	-6.9	48,000	1229	475	447	-28.7	62,300
1115	1385	762	-6.4	40,400	1230	466	1282	-29.0	20,400 14,400
1116	1078	816	-10.6	38,000	1231 1232	759 1324	1461 1170	-17.4 -7.2	24,200
1117	975	787 033	-12.6 -8.7	39,300 33,100	1232	1583	1005	-7.2	30,300
1118	1202 1022	933 1076	-8.7 -11.6	27,600	1234	1865	809	-0.6	38,200
1119 1120	1905	616	-0.3	48,300	1235	1812	817	-1.0	37,900
1121	1512	1301	-4.5	19,700	1236	1411	703	-6.0	43,400
1122	1114	677	-9.9	44,700	1237	1392	682	-6.3	44,500
1123	1464	452	-5.1	61,700	1238	794	410	-16.4	66,900
1125	1048	857	-11.1	36,200	1239	769	407	-17.1	67,300
1126	1122	802	-9.8	38,600	1240	740	406	-17.9	67,500
1128	1722	892	-2.1	34,700	1241	743	511	-17.8	55,900 56,000
1133	1098	825	-10.2	37,500	1242	713	510 509	-18.7 -19.6	56,000 56,100
1139	1830	569	-0.8	51,400 23,800	1243 1244	682 663	504	-19.6	56,500
1147 1148	764 1968	1182 724	-17.3 >0.0	23,800 42,300	1244	565	582	-24.4	50,500
1140	1500	124	20.0	-2,000					

MSN X Y CPKpl SDSMW 1246 547 577 -25.3 50,800 1247 530 576 -26.3 50,900 1249 516 572 -27.0 51,200 1250 973 536 -12.7 53,900 1251 607 532 -22.4 54,200 1252 665 529 -20.2 54,400 1253 899 766 -14.1 40,200 1254 1311 746 -7.4 41,200 1255 1300 761 -7.5 40,400 1257 1938 712 0.0 42,600 1258 1806 718 -1.0 42,600 1259 1727 715 -2.0 42,700 1260 1629 713 -3.0 42,600 1261 1555 717 -4.0 42,600 1262 1468 717 -5.0					
1247 530 576 -26.3 50,900 1249 516 572 -27.0 51,200 1250 973 536 -12.7 53,900 1251 607 532 -22.4 54,200 1252 665 529 -20.2 54,400 1253 899 766 -14.1 40,200 1255 1300 761 -7.5 40,400 1257 1938 712 0.0 42,900 1258 1806 718 -1.0 42,600 1259 1727 715 -2.0 42,700 1260 1629 713 -3.0 42,600 1261 1555 717 -4.0 42,600 1262 1468 717 -5.0 42,600 1263 1413 722 -6.0 42,400 1264 1340 717 -7.0 42,600 1266 1182 720 -9.0	MSN	x	Y	CPKpI	SDSMW
1249 516 572 -27.0 51,200 1250 973 536 -12.7 53,900 1251 607 532 -22.4 54,200 1252 665 529 -20.2 54,400 1253 899 766 -14.1 40,200 1255 1300 761 -7.5 40,400 1255 1300 761 -7.5 40,400 1257 1938 712 0.0 42,600 1258 1806 718 -1.0 42,600 1258 1806 718 -1.0 42,600 1259 1727 715 -2.0 42,700 1260 1629 713 -3.0 42,600 1261 1555 717 -4.0 42,600 1262 1468 717 -5.0 42,600 1263 1413 722 -6.0 42,600 1264 1340 717 -10.0 <td>1246</td> <td>547</td> <td>577</td> <td>-25.3</td> <td>50,800</td>	1246	547	577	-25.3	50,800
1250 973 536 -12.7 53,900 1251 607 532 -22.4 54,200 1252 665 529 -20.2 54,400 1253 899 766 -14.1 40,200 1254 1311 746 -7.4 41,200 1255 1300 761 -7.5 40,400 1257 1938 712 0.0 42,900 1258 1806 718 -1.0 42,600 1259 1727 715 -2.0 42,700 1260 1629 713 -3.0 42,600 1261 1555 717 -4.0 42,600 1262 1468 717 -5.0 42,600 1264 1340 717 -7.0 42,600 1265 1263 717 -8.0 42,600 1266 1182 720 -9.0 42,500 1267 1110 717 -10.0 <td>1247</td> <td>530</td> <td>576</td> <td>-26.3</td> <td>50,900</td>	1247	530	576	-26.3	50,900
1251 607 532 -22.4 54,200 1252 665 529 -20.2 54,400 1253 899 766 -14.1 40,200 1254 1311 746 -7.4 41,200 1255 1300 761 -7.5 40,400 1257 1938 712 0.0 42,900 1258 1806 718 -1.0 42,600 1259 1727 715 -2.0 42,700 1260 1629 713 -3.0 42,600 1261 1555 717 -4.0 42,600 1262 1468 717 -5.0 42,600 1264 1340 717 -7.0 42,600 1265 1263 717 -8.0 42,600 1266 1182 720 -9.0 42,600 1268 1055 717 -11.0 42,600 1269 110 717 -10.0 <td>1249</td> <td>516</td> <td>572</td> <td>-27.0</td> <td>51,200</td>	1249	516	572	-27.0	51,200
1252 665 529 -20.2 54,400 1253 899 766 -14.1 40,200 1254 1311 746 -7.4 41,200 1255 1300 761 -7.5 40,400 1257 1938 712 0.0 42,600 1258 1806 718 -1.0 42,600 1259 1727 715 -2.0 42,700 1260 1629 713 -3.0 42,600 1261 1555 717 -4.0 42,600 1262 1468 717 -5.0 42,600 1263 1413 722 -6.0 42,400 1265 1263 717 -8.0 42,600 1266 1182 720 -9.0 42,500 1268 1055 717 -11.0 42,600 1270 959 715 -13.0 42,700 1271 905 712 -14.0 <td>1250</td> <td>973</td> <td>536</td> <td>-12.7</td> <td>53,900</td>	1250	973	536	-12.7	53,900
1253 899 766 -14.1 40,200 1254 1311 746 -7.4 41,200 1255 1300 761 -7.5 40,400 1257 1938 712 0.0 42,900 1258 1806 718 -1.0 42,600 1259 1727 715 -2.0 42,700 1260 1629 713 -3.0 42,800 1261 1555 717 -4.0 42,600 1262 1468 717 -5.0 42,600 1263 1413 722 -6.0 42,400 1264 1340 717 -7.0 42,600 1265 1263 717 -8.0 42,600 1266 1182 720 -9.0 42,500 1267 1110 717 -10.0 42,600 1268 1055 717 -11.0 42,600 1267 1110 717 -10.0 </td <td>1251</td> <td>607</td> <td>532</td> <td>-22.4</td> <td>54,200</td>	1251	607	532	-22.4	54,200
1254 1311 746 -7.4 41,200 1255 1300 761 -7.5 40,400 1257 1938 712 0.0 42,900 1258 1806 718 -1.0 42,600 1259 1727 715 -2.0 42,700 1260 1629 713 -3.0 42,600 1261 1555 717 -4.0 42,600 1262 1468 717 -5.0 42,600 1263 1413 722 -6.0 42,600 1264 1340 717 -7.0 42,600 1265 1263 717 -8.0 42,600 1266 1182 720 -9.0 42,500 1267 1110 717 -10.0 42,600 1268 1055 717 -11.0 42,600 1269 999 717 -12.0 42,600 1271 905 712 -14.0 <td>1252</td> <td>665</td> <td>529</td> <td>-20.2</td> <td>54,400</td>	1252	665	529	-20.2	54,400
1255 1300 761 -7.5 40,400 1257 1938 712 0.0 42,900 1258 1806 718 -1.0 42,600 1259 1727 715 -2.0 42,700 1260 1629 713 -3.0 42,800 1261 1555 717 -4.0 42,600 1262 1468 717 -5.0 42,600 1263 1413 722 -6.0 42,600 1264 1340 717 -7.0 42,600 1265 1263 717 -8.0 42,600 1265 1263 717 -10.0 42,600 1266 1182 720 -9.0 42,600 1267 1110 717 -10.0 42,600 1268 1055 717 -11.0 42,600 1270 959 715 -13.0 42,700 1271 905 712 -14.0 </td <td>1253</td> <td>899</td> <td>766</td> <td>-14.1</td> <td>40,200</td>	1253	899	766	-14.1	40,200
1257 1938 712 0.0 42,900 1258 1806 718 -1.0 42,600 1259 1727 715 -2.0 42,700 1260 1629 713 -3.0 42,800 1261 1555 717 -4.0 42,600 1262 1468 717 -5.0 42,600 1263 1413 722 -6.0 42,400 1264 1340 717 -7.0 42,600 1265 1263 717 -8.0 42,600 1265 1263 717 -10.0 42,500 1266 1182 720 -9.0 42,500 1267 1110 717 -10.0 42,600 1268 1055 717 -11.0 42,600 1270 959 715 -13.0 42,700 1271 905 712 -14.0 42,900 1271 905 712 -14.0 </td <td>1254</td> <td>1311</td> <td>746</td> <td>-7.4</td> <td>41,200</td>	1254	1311	746	-7.4	41,200
1257 1938 712 0.0 42,900 1258 1806 718 -1.0 42,600 1259 1727 715 -2.0 42,700 1260 1629 713 -3.0 42,800 1261 1555 717 -4.0 42,600 1262 1468 717 -5.0 42,600 1263 1413 722 -6.0 42,400 1264 1340 717 -7.0 42,600 1265 1263 717 -8.0 42,600 1266 1182 720 -9.0 42,500 1267 1110 717 -10.0 42,600 1268 1055 717 -11.0 42,600 1269 999 717 -12.0 42,600 1270 959 715 -13.0 42,700 1271 905 712 -14.0 42,900 1271 957 714 -15.0 <td></td> <td>1300</td> <td>761</td> <td>-7.5</td> <td>40,400</td>		1300	761	-7.5	40,400
1259 1727 715 -2.0 42,700 1260 1629 713 -3.0 42,800 1261 1555 717 -4.0 42,600 1262 1468 717 -5.0 42,600 1263 1413 722 -6.0 42,400 1265 1263 717 -8.0 42,600 1266 1182 720 -9.0 42,600 1267 1110 717 -10.0 42,600 1268 1055 717 -11.0 42,600 1269 999 717 -12.0 42,600 1270 959 715 -13.0 42,700 1271 905 712 -14.0 42,900 1272 857 714 -15.0 42,800 1273 810 705 -16.0 43,300 1274 774 711 -17.0 42,900 1277 737 708 -18.0 </td <td>1257</td> <td>1938</td> <td>712</td> <td>0.0</td> <td></td>	1257	1938	712	0.0	
1259 1727 715 -2.0 42,700 1260 1629 713 -3.0 42,800 1261 1555 717 -4.0 42,600 1262 1468 717 -5.0 42,600 1263 1413 722 -6.0 42,400 1265 1263 717 -8.0 42,600 1266 1182 720 -9.0 42,500 1267 1110 717 -10.0 42,600 1268 1055 717 -11.0 42,600 1268 1055 717 -11.0 42,600 1268 1055 717 -11.0 42,600 1270 959 715 -13.0 42,700 1271 905 712 -14.0 42,900 1271 905 712 -14.0 42,900 1272 857 714 -15.0 42,900 1273 810 705 -16.0		1806		-1.0	42,600
1260 1629 713 -3.0 42,800 1261 1555 717 -4.0 42,600 1262 1468 717 -5.0 42,600 1263 1413 722 -6.0 42,400 1264 1340 717 -7.0 42,600 1265 1263 717 -8.0 42,600 1266 1182 720 -9.0 42,500 1267 1110 717 -10.0 42,600 1268 1055 717 -11.0 42,600 1269 999 717 -12.0 42,600 1270 959 715 -13.0 42,700 1271 905 712 -14.0 42,900 1272 857 714 -15.0 42,800 1273 810 705 -16.0 43,300 1274 774 711 -17.0 42,900 1278 702 711 -19.0 </td <td></td> <td>1727</td> <td>715</td> <td>-2.0</td> <td>42,700</td>		1727	715	-2.0	42,700
1261 1555 717 -4.0 42,600 1262 1468 717 -5.0 42,600 1263 1413 722 -6.0 42,400 1264 1340 717 -7.0 42,600 1265 1263 717 -8.0 42,600 1266 1182 720 -9.0 42,500 1267 1110 717 -10.0 42,600 1268 1055 717 -11.0 42,600 1269 999 717 -12.0 42,600 1270 959 715 -13.0 42,700 1271 905 712 -14.0 42,900 1271 905 712 -14.0 42,900 1272 857 714 -15.0 42,800 1273 810 705 -16.0 43,300 1274 774 711 -17.0 42,900 1273 707 708 -18.0 </td <td></td> <td></td> <td>713</td> <td>-3.0</td> <td>42,800</td>			713	-3.0	42,800
1262 1468 717 -5.0 42,600 1263 1413 722 -6.0 42,400 1264 1340 717 -7.0 42,600 1265 1263 717 -8.0 42,600 1266 1182 720 -9.0 42,600 1268 1055 717 -11.0 42,600 1269 999 717 -12.0 42,600 1270 959 715 -13.0 42,700 1271 905 712 -14.0 42,900 1271 905 712 -14.0 42,900 1272 857 714 -15.0 42,800 1273 810 705 -16.0 43,300 1274 774 711 -17.0 42,900 1278 702 711 -19.0 42,900 1279 671 710 -20.0 43,000 1281 617 707 -22.0 <td></td> <td></td> <td>717</td> <td>-4.0</td> <td>42,600</td>			717	-4.0	42,600
1263 1413 722 -6.0 42,400 1264 1340 717 -7.0 42,600 1265 1263 717 -8.0 42,600 1266 1182 720 -9.0 42,500 1267 1110 717 -10.0 42,600 1268 1055 717 -11.0 42,600 1270 959 715 -13.0 42,700 1271 905 712 -14.0 42,800 1272 857 714 -15.0 42,800 1273 810 705 -16.0 43,300 1274 774 711 -17.0 42,900 1277 737 708 -18.0 43,100 1278 702 711 -19.0 42,900 1279 671 710 -20.0 43,000 1280 645 710 -21.0 43,000 1281 617 707 -22.0 </td <td></td> <td></td> <td>717</td> <td>-5.0</td> <td>42,600</td>			717	-5.0	42,600
1264 1340 717 -7.0 42,600 1265 1263 717 -8.0 42,600 1266 1182 720 -9.0 42,500 1267 1110 717 -10.0 42,600 1268 1055 717 -11.0 42,600 1269 999 717 -12.0 42,600 1270 959 715 -13.0 42,700 1271 905 712 -14.0 42,900 1272 857 714 -15.0 42,800 1273 810 705 -16.0 43,300 1274 774 711 -17.0 42,900 1277 737 708 -18.0 43,100 1278 702 711 -19.0 42,900 1279 671 710 -20.0 43,000 1280 645 710 -21.0 43,000 1281 617 707 -22.0 </td <td></td> <td>1413</td> <td>722</td> <td>-6.0</td> <td>42,400</td>		1413	722	-6.0	42,400
1265 1263 717 -8.0 42,600 1266 1182 720 -9.0 42,500 1267 1110 717 -10.0 42,600 1268 1055 717 -11.0 42,600 1269 999 717 -12.0 42,600 1270 959 715 -13.0 42,700 1271 905 712 -14.0 42,900 1272 857 714 -15.0 42,800 1273 810 705 -16.0 43,300 1274 774 711 -17.0 42,900 1277 737 708 -18.0 43,100 1278 702 711 -19.0 42,900 1279 671 710 -20.0 43,000 1280 645 710 -21.0 43,000 1281 617 707 -22.0 43,300 1282 595 704 -23.0 </td <td></td> <td>1340</td> <td>717</td> <td>-7.0</td> <td>42,600</td>		1340	717	-7.0	42,600
1266 1182 720 -9.0 42,500 1267 1110 717 -10.0 42,600 1268 1055 717 -11.0 42,600 1269 999 717 -12.0 42,600 1270 959 715 -13.0 42,700 1271 905 712 -14.0 42,900 1272 857 714 -15.0 42,800 1273 810 705 -16.0 43,300 1274 774 711 -17.0 42,900 1277 737 708 -18.0 43,100 1278 702 711 -19.0 42,900 1279 671 710 -20.0 43,000 1280 645 710 -21.0 43,000 1281 617 707 -22.0 43,100 1282 595 704 -23.0 43,300 1284 552 695 -25.0 </td <td></td> <td></td> <td>717</td> <td>-8.0</td> <td>42,600</td>			717	-8.0	42,600
1267 1110 717 -10.0 42,600 1268 1055 717 -11.0 42,600 1269 999 717 -12.0 42,600 1270 959 715 -13.0 42,700 1271 905 712 -14.0 42,900 1272 857 714 -15.0 42,800 1273 810 705 -16.0 43,300 1274 774 711 -17.0 42,900 1277 737 708 -18.0 43,100 1278 702 711 -19.0 42,900 1279 671 710 -20.0 43,000 1280 645 710 -21.0 43,000 1281 617 707 -22.0 43,100 1282 595 704 -23.0 43,300 1283 573 700 -24.0 43,500 1284 552 695 -25.0 </td <td></td> <td>1182</td> <td>720</td> <td>-9.0</td> <td>42,500</td>		1182	720	-9.0	42,500
1269 999 717 -12.0 42,600 1270 959 715 -13.0 42,700 1271 905 712 -14.0 42,900 1272 857 714 -15.0 42,800 1273 810 705 -16.0 43,300 1274 774 711 -17.0 42,900 1278 702 711 -19.0 42,900 1279 671 710 -20.0 43,000 1280 645 710 -21.0 43,000 1281 617 707 -22.0 43,100 1282 595 704 -23.0 43,300 1283 573 700 -24.0 43,500 1284 552 695 -25.0 43,700 1285 536 694 -26.0 43,800 1286 515 687 -27.0 44,200 1287 496 683 -28.0 <td></td> <td>1110</td> <td>717</td> <td>-10.0</td> <td>42,600</td>		1110	717	-10.0	42,600
1270 959 715 -13.0 42,700 1271 905 712 -14.0 42,900 1272 857 714 -15.0 42,800 1273 810 705 -16.0 43,300 1274 774 711 -17.0 42,900 1277 737 708 -18.0 43,100 1278 702 711 -19.0 42,900 1279 671 710 -20.0 43,000 1280 645 710 -21.0 43,000 1281 617 707 -22.0 43,100 1282 595 704 -23.0 43,300 1283 573 700 -24.0 43,500 1284 552 695 -25.0 43,700 1285 536 694 -26.0 43,800 1286 515 687 -27.0 44,200 1287 496 683 -28.0 <td>1268</td> <td>1055</td> <td>717</td> <td>-11.0</td> <td>42,600</td>	1268	1055	717	-11.0	42,600
1271 905 712 -14.0 42,900 1272 857 714 -15.0 42,800 1273 810 705 -16.0 43,300 1274 774 711 -17.0 42,900 1277 737 708 -18.0 43,100 1278 702 711 -19.0 42,900 1279 671 710 -20.0 43,000 1280 645 710 -21.0 43,000 1281 617 707 -22.0 43,100 1282 595 704 -23.0 43,300 1283 573 700 -24.0 43,500 1284 552 695 -25.0 43,700 1285 536 694 -26.0 43,800 1286 515 687 -27.0 44,200 1287 496 683 -28.0 44,400 1288 467 669 -29.0 <td>1269</td> <td>999</td> <td>717</td> <td>-12.0</td> <td>42,600</td>	1269	999	717	-12.0	42,600
1272 857 714 -15.0 42,800 1273 810 705 -16.0 43,300 1274 774 711 -17.0 42,900 1277 737 708 -18.0 43,100 1278 702 711 -19.0 42,900 1279 671 710 -20.0 43,000 1280 645 710 -21.0 43,000 1281 617 707 -22.0 43,100 1282 595 704 -23.0 43,300 1283 573 700 -24.0 43,500 1284 552 695 -25.0 43,700 1284 552 695 -25.0 43,800 1285 536 694 -26.0 43,800 1287 496 683 -28.0 44,400 1288 467 669 -29.0 45,300 1290 427 655 -31.0 <td>1270</td> <td>959</td> <td>715</td> <td>-13.0</td> <td>42,700</td>	1270	959	715	-13.0	42,700
1273 810 705 -16.0 43,300 1274 774 711 -17.0 42,900 1277 737 708 -18.0 43,100 1278 702 711 -19.0 42,900 1279 671 710 -20.0 43,000 1280 645 710 -21.0 43,000 1281 617 707 -22.0 43,100 1282 595 704 -23.0 43,300 1283 573 700 -24.0 43,500 1284 552 695 -25.0 43,700 1284 552 695 -25.0 43,800 1285 536 694 -26.0 43,800 1286 515 687 -27.0 44,200 1287 496 683 -28.0 44,400 1288 467 669 -29.0 45,300 1290 427 655 -31.0 <td>1271</td> <td>905</td> <td>712</td> <td>-14.0</td> <td>42,900</td>	1271	905	712	-14.0	42,900
1274 774 711 -17.0 42,900 1277 737 708 -18.0 43,100 1278 702 711 -19.0 42,900 1279 671 710 -20.0 43,000 1280 645 710 -21.0 43,000 1281 617 707 -22.0 43,100 1282 595 704 -23.0 43,300 1283 573 700 -24.0 43,500 1284 552 695 -25.0 43,700 1285 536 694 -26.0 43,800 1286 515 687 -27.0 44,200 1287 496 683 -28.0 44,200 1288 467 669 -29.0 45,200 1289 447 667 -30.9 45,300 1291 412 655 -31.0 45,900 1292 397 652 -33.0 <td>1272</td> <td>857</td> <td>714</td> <td>-15.0</td> <td>42,800</td>	1272	857	714	-15.0	42,800
1277 737 708 -18.0 43,100 1278 702 711 -19.0 42,900 1279 671 710 -20.0 43,000 1280 645 710 -21.0 43,000 1281 617 707 -22.0 43,100 1282 595 704 -23.0 43,500 1283 573 700 -24.0 43,500 1284 552 695 -25.0 43,700 1285 536 694 -26.0 43,800 1286 515 687 -27.0 44,200 1287 496 683 -28.0 44,400 1288 467 669 -29.0 45,200 1289 447 667 -30.9 45,300 1291 412 655 -31.0 45,900 1292 397 652 -33.0 46,100 1293 381 654 -34.0 <td>1273</td> <td>810</td> <td>705</td> <td>-16.0</td> <td>43,300</td>	1273	810	705	-16.0	43,300
1278 702 711 -19.0 42,900 1279 671 710 -20.0 43,000 1280 645 710 -21.0 43,000 1281 617 707 -22.0 43,100 1282 595 704 -23.0 43,300 1283 573 700 -24.0 43,500 1284 552 695 -25.0 43,700 1285 536 694 -26.0 43,800 1286 515 687 -27.0 44,200 1287 496 683 -28.0 44,400 1288 467 669 -29.0 45,200 1289 447 667 -30.9 45,300 1290 427 655 -31.0 45,900 1291 412 655 -32.0 46,100 1293 381 654 -34.0 46,000 1294 365 653 -35.0 <td>1274</td> <td>774</td> <td>711</td> <td>-17.0</td> <td>42,900</td>	1274	774	711	-17.0	42,900
1279 671 710 -20.0 43,000 1280 645 710 -21.0 43,000 1281 617 707 -22.0 43,100 1282 595 704 -23.0 43,300 1283 573 700 -24.0 43,500 1284 552 695 -25.0 43,700 1285 536 694 -26.0 43,800 1286 515 687 -27.0 44,200 1287 496 683 -28.0 44,400 1288 467 669 -29.0 45,200 1289 447 667 -30.9 45,300 1290 427 655 -31.0 45,900 1291 412 655 -32.0 45,900 1293 381 654 -34.0 46,100 1294 365 653 -35.0 46,100	1277	737	708	-18.0	43,100
1280 645 710 -21.0 43,000 1281 617 707 -22.0 43,100 1282 595 704 -23.0 43,300 1283 573 700 -24.0 43,500 1284 552 695 -25.0 43,700 1285 536 694 -26.0 43,800 1286 515 687 -27.0 44,200 1287 496 683 -28.0 44,400 1288 467 669 -29.0 45,200 1289 447 667 -30.9 45,300 1290 427 655 -31.0 45,900 1291 412 655 -32.0 45,900 1292 397 652 -33.0 46,100 1293 381 654 -34.0 46,000 1294 365 653 -35.0 46,100	1278	702	711	-19.0	42,900
1281 617 707 -22.0 43,100 1282 595 704 -23.0 43,300 1283 573 700 -24.0 43,500 1284 552 695 -25.0 43,700 1285 536 694 -26.0 43,800 1286 515 687 -27.0 44,200 1287 496 683 -28.0 44,400 1288 467 669 -29.0 45,200 1289 447 667 -30.9 45,300 1290 427 655 -31.0 45,900 1291 412 655 -32.0 45,900 1292 397 652 -33.0 46,100 1293 381 654 -34.0 46,000 1294 365 653 -35.0 46,100	1279	671	710	-20.0	43,000
1282 595 704 -23.0 43,300 1283 573 700 -24.0 43,500 1284 552 695 -25.0 43,700 1285 536 694 -26.0 43,800 1286 515 687 -27.0 44,200 1287 496 683 -28.0 44,400 1288 467 669 -29.0 45,200 1289 447 667 -30.9 45,300 1290 427 655 -31.0 45,900 1291 412 655 -32.0 45,900 1292 397 652 -33.0 46,100 1293 381 654 -34.0 46,000 1294 365 653 -35.0 46,100	1280	645	710	-21.0	43,000
1283 573 700 -24.0 43,500 1284 552 695 -25.0 43,700 1285 536 694 -26.0 43,800 1286 515 687 -27.0 44,200 1287 496 683 -28.0 44,400 1288 467 669 -29.0 45,200 1289 447 667 -30.9 45,300 1290 427 655 -31.0 45,900 1291 412 655 -32.0 45,900 1292 397 652 -33.0 46,100 1293 381 654 -34.0 46,000 1294 365 653 -35.0 46,100	1281	617	707	-22.0	43,100
1284 552 695 -25.0 43,700 1285 536 694 -26.0 43,800 1286 515 687 -27.0 44,200 1287 496 683 -28.0 44,400 1288 467 669 -29.0 45,200 1289 447 667 -30.9 45,300 1290 427 655 -31.0 45,900 1291 412 655 -32.0 45,900 1292 397 652 -33.0 46,100 1293 381 654 -34.0 46,000 1294 365 653 -35.0 46,100	1282	595	704	-23.0	43,300
1285 536 694 -26.0 43,800 1286 515 687 -27.0 44,200 1287 496 683 -28.0 44,400 1288 467 669 -29.0 45,200 1299 447 667 -30.9 45,300 1290 427 655 -31.0 45,900 1291 412 655 -32.0 45,900 1292 397 652 -33.0 46,100 1293 381 654 -34.0 46,000 1294 365 653 -35.0 46,100	1283	573	700	-24.0	43,500
1286 515 687 -27.0 44,200 1287 496 683 -28.0 44,400 1288 467 669 -29.0 45,200 1289 447 667 -30.9 45,300 1290 427 655 -31.0 45,900 1291 412 655 -32.0 45,900 1292 397 652 -33.0 46,100 1293 381 654 -34.0 46,000 1294 365 653 -35.0 46,100	1284	552	695	-25.0	43,700
1287 496 683 -28.0 44,400 1288 467 669 -29.0 45,200 1289 447 667 -30.9 45,300 1290 427 655 -31.0 45,900 1291 412 655 -32.0 45,900 1292 397 652 -33.0 46,100 1293 381 654 -34.0 46,000 1294 365 653 -35.0 46,100	1285	536	694	-26.0	43,800
1288 467 669 -29.0 45,200 1289 447 667 -30.9 45,300 1290 427 655 -31.0 45,900 1291 412 655 -32.0 45,900 1292 397 652 -33.0 46,100 1293 381 654 -34.0 46,000 1294 365 653 -35.0 46,100	1286	515	687	-27.0	44,200
1289 447 667 -30.9 45,300 1290 427 655 -31.0 45,900 1291 412 655 -32.0 45,900 1292 397 652 -33.0 46,100 1293 381 654 -34.0 46,000 1294 365 653 -35.0 46,100	1287	496	683		44,400
1290 427 655 -31.0 45,900 1291 412 655 -32.0 45,900 1292 397 652 -33.0 46,100 1293 381 654 -34.0 46,000 1294 365 653 -35.0 46,100	1288	467	669		
1291 412 655 -32.0 45,900 1292 397 652 -33.0 46,100 1293 381 654 -34.0 46,000 1294 365 653 -35.0 46,100					•
1292 397 652 -33.0 46,100 1293 381 654 -34.0 46,000 1294 365 653 -35.0 46,100					
1293 381 654 -34.0 46,000 1294 365 653 -35.0 46,100					
1294 365 653 -35.0 46,100					
. =					
1295 348 653 <-35.0 46,100					
	1295	348	653	<-35.0	46,100
			-		

POP name	Protein name	MSN's	Basis for identification
IDS:3_ALPHA_HDDH	3-a-hydroxysteroid-dihydrodiol- dehydrogenase, an enzyme of	137, 159	Pure protein and antibody provided by Dr. T.M. Penning, Department of Pharmacology, School
IDS:ACTIN_BETA	steroid metabolism β cellular actin, a cytoskeletal protein	38	of Medicine, University of Pennsylvania. Homologous position with respect to other mammalian
IDS:ACTIN_GAMMA	γ cellular actin, a cytoskeletal protein	89	systems Homologous position with respect to other mammalian
IDS:ALBUMIN IDS:APO_A-I	Serum albumin, mature form. Apo A-I plasma lipoprotein, mature form	21, 28, 33 236, 463	Predominance in rat plasma Presence in rat plasma, regulation by some lipid-
IDS:CALMODULIN	(tentative). Calmodulin, an acidic cytosolic calcium-	123, 649	towering arings Homologous position with respect to other mammalian
IDS:CATALASE	binding protein Catalase (peroxisomal)	54, 61, 106	systems Presence in purified peroxisomes, similarity in position
IDS:CPKSPOTS	Spots contributed by the CPK charge	1257 - 1295	
IDS:CPS	standards (not rat liver proteins) Carbamoyl phosphate synthase	114, 157, 167, 174, 1184, 1185, 1186, 1222	Pure protein provided by Dr. Margaret Marshall, Department of Pharmacology, Medical School,
IDS:CYTOCHROME_B5	Cytochrame b5	87, 477	University of Wisconsin - Madison. Pure protein provided by Dr. Andrew Parkinson, Department of Pharmacology, Toxicology and Therapeutics, University of Kansas Medical
IDS:FABP-L	Liver fatty-acid binding protein	227	Center Pure protein provided by Dr. Nathan Bass, Department of Medicine, University of California School of
IDS;HMG-COA_SYNTHASE	Cytosolic HMG-CoA Symhase	133, 144, 235, 413	Antibody provided by Dr. Michael Greenspan, Merck Sharp & Dohme Research Laboratories,
IDS:LAMIN_B	Lamin B, a nuclear protein	415, 734	Homologous position with respect to other mammalian
IDS:MITCON:1	Mitcon:1 (F1 ATPase B subunit), a	17, 49, 71, 340, 1245, 1246, 1247, 1249	Homologous position with respect to other mammalian
IDS:MITCON:2	Mitcon:2, a mitochondrial matrix stress	15, 25, 110, 1241, 1242, 1243, 1244	Homologous position with respect to other mammallan
IDS:MITCON:3	protein equivalent to E. Mitcon:3, a mitcohondrial matrix stress	18, 35, 226, 600, 1238, 1239, 1240	Homologous position with respect to other mammalian
IDS:NADPH_P450_RED	protein, likely analog of NADPH cytochrome P-450 reductase, frequently co-induced with P-450's	175, 251, 812	Systems, presence in minociprorial Pure profein provided by Dr. Andrew Parkinson, Department of Pharmacology, Toxicology and Therapeutics, University of Kansas Medical
IDS:PDI	Protein disulphide isomerase 1	168, 1170, 1171, 1172	Sequence information obtained by R.M. Van Frank,
IDS:PLASMA_PROTEINS	Rat plasma proteins observed in liver	21, 28, 33, 44, 72, 102, 115, 197, 236, 246, 248, 257, 293, 332, 347, 364, 369, 419, 432, 463, 463, 463, 562, 665, 623, 666, 667, 725, 738, 738, 738, 738, 738, 738, 738, 738	Plasma coelectrophoresis studies
IDS:PRO-ALBUMIN	Serum albumin precursor	736, 730, 803, 303, 320 47, 93	Relative position to mature albumin, presence in micro-
IDS:PYRCARBOX IDS:SOD	Pyruvate carboxylase Superoxide dismutase	179, 1180, 1181, 1182, 1183 135	Paviica, R.J., et al., BBA (1990) 1022 115-125. Sequence information obtained by R.M. Van Frank,
IDS:TUBULIN_ALPHA	lpha tubulin, a cytoskeletal protein	56, 132, 1224, 1252	Homologous position with respect to other mammalian
IDS-TUBBILIN BETA	R tubulin a Adoskolotal protein	50, 1225, 1226, 1251	Homologous position with respect to other mammalian


Table 3. Computed pl's of two sets of carbamylated protein standards: Rabbit muscle CPK and human hemoglobin (Hb)

	Protein Name	PIR Name	#ASP 3.9	#GLU 4.1	#HIS 6.0	#LYS 10.8	#ARG 12.5	NH2- 7.0	Calc	Real CPK
0	Rabbit muscle CPK	KIRBCM	28	27	17	34	18	1	6.84	0.0
-1			28	27	17	33	18	1	6.67	-1
-2			28	27	17	32	18	1	6.54	-2
-3			28	27	17	31	18	1	6.42	-3
-4			28	27	17	30	18	1	6.31	-4
-5			28	27	17	29	18	1	6.21	-5
-6			28	27	17	28	18	1	6.12	-6
-7			28	27	17	27	18	1	6.03	-7 -8
-8			28	27	17	26	18	1	5.94	-8 -9
-9 10			28 28	27 27	17 17	25 24	18 18	1	5.85 5.76	-10
-10 -11			28	27	17	23	18	1	5.67	-11
-12			28	27	17	22	18	•	5.58	-12
-13			28	27	17	21	18	i	5.48	-13
-14			28	27	17	20	18	1	5.39	-14
-15			28	27	17	19	18	1	5.29	-15
-16			28	27	17	18	18	1	5.20	-16
-17			28	27	17	17	18	1	5.12	-17
-18			28	27	17	16	18	1	5.04	-18
-19			28	27	17	15	18	1	4.96	-19
-20			28	27	17	14	18	1	4.89	-20
-21			28	27	17	13	18	1	4.83	-21
-22			28	27	17	12	18	1	4.77	-22
-23			28	27	17	11	18	1	4.71	-23
-24			28	27	17	10	18	1	4.66	-24
-25			28	27	17	9	18	1	4.61	-25
-26			28 28	27 27	17 17	8 7	18 18	1	4.56 4.52	-26 -27
-27 -28			28	27	17	6	18	1	4.48	-28
-29			28	27	17	5	18	1	4.44	-29
-30			28	27	17	4	18	1	4.40	-30
-31			28	27	17	3	18	1	4.36	-31
-32			28	27	17	2	18	1	4.32	-32
-33			28	27	17	1	18	1	4.29	-33
-34			28	27	17	0	18	1	4.25	-34
-35			28	27	17	0	18	0	4.22	-35
0	Hb-beta, human	HBHU	7	8	9	11	3	1	7.18	
-1	·		7	8	9	10	3	1	6.79	
-2			.7	8	9	9	3	1	6.53	-1.8
-3			7 7	8	9	8	3	1	6.32	-3.2
-4			7	8	9	7	3	1	6.13	-5.3
-5	•		7	8	9	6	3	1	5.96	-7.2
-6			7	8	9	5	3	1	5.78	-10.0
-7 .			7	8	9	4	3	1	5.59	-12.3
-8			7	8	9	3	3	1	5.37	-15.5
-9			7	8	9	2	3	1	5.14	-18.0
-10			7	8	9	1	3	1	4.91	-21.0
-11			7	8	9	0	3	1	4.71	-25.5
-12			7	8	9	0	3	0	4.54	-27.2

Table 4. Computed pl's of some known proteins related to measured CPK pl's


audio 4.		PIR Name	#ASP 3.9	#GLU 4.1	#HIS 6.0		#ARG 12.5	Calc pl	Real CPK
	Protein Name	MDDOM	28	27	17	34	18	6.84	0.0
0	Creatine phospho kinase (CPK), rabbit muscle	KIRBCM FZRTL	26 5	13	2	16	2	7.83	-3.0
1	Fatty acid-binding protein, rat nepatic	MGHUB2		8	4	8	5	6.09	-5.0
2	h2-microglobulin, human	SYRTCA	72	96	28	_	56	5.97	-5.5
3	Carbamovi-phosphate synthase, rat	ABRTS	32	57	15		27	5.98	-6.2
4	Proalbumin (serum albumin precursor), rat	ABRTS	32		15		24	5.71	-9.0
5	Sozum albumin, rat	A26810	8	-	10		4	5.91	-9.2
6	Superoxid dismutase (Cu-Zn, SOD), rat	A28807	34		9		21	5.92	-9.2
7	Phospholipase C, phophoinositide-specific (?), rat	ABHUS	36		16			5.70	-11.9
8	Albumin, human	A24700	18	_	. 6			5.32	-13.7
9	Apo A-I lipoprotein, rat	LPHUA1	16		6		17	5.35	-14.3
10	proApo A-I lipoprotein, human	RDRTO4			_		36	5.07	-15.6
11	NADPH cytochrome P-450 reductase, rat	VAHU	18					5.04	-16.9
12	Retinol binding protein, human	ATRTC	23				18	5.06	-17.2
13	Actin beta, rat	ATRIC	20				18	5.07	-16.8
14	Actin gamma, rat	LPHUA1	16					5.10	-17.5
15	Apo A-I lipoprotein, human	LPHUA4	20				24	4.88	-19.7
16	Apo A-IV lipoprotein, human	UBRTA	27					4.66	-19.8
17	Tubulin alpha, rat	PWBOB	25					4.80	-21.0
18	F1ATPase beta, bovine	UBPGB	26						-22.5
19	Tubulin beta, pig	ISRTSS	43						-25.0
20	Protein disulphide isomerase (PDI), rat hepatic	CBRT5	10						-26.0
21	Cytochrome b5, rat	LPHUC2						4.44	-30.5
22	Apo C-II lipoprotein, human	LFHUUZ		, ,					
	Amino acid pl assumed in calulation:		3.9	4.1	6.0	10.8	12.5		

An International Journal

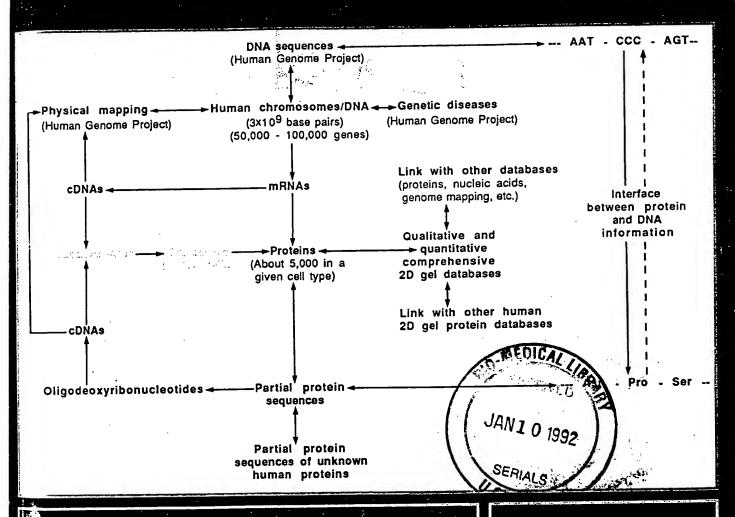

TWO-DIMENSIONAL GEL PROTEIN DATABASES

Editor: J. E. Celis

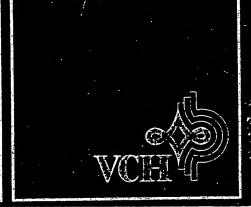

ELECTROPHORESIS

Indexed in: BIOSI Current Contents, MEDLAR ISSN 0173-083 ELCTDN 12 (11) 763–996 (1991

An International Journal

TWO-DIMENSIONAL GEL PROTEIN DATABASES Editor: J. E. Celis

Editorial

The master two-dimensional gel database of human AMA cell proteins: Towards J. E. Celis, H. Leffers, linking protein and genome sequence and mapping information (Update 1991) H. H. Rasmussen, P. Madsen, B. Honoré, B. Gesser, K. Dejgaard, E. Olsen, G. P. Ratz, J. B. Lauridsen, B. Basse, A. H. Andersen, E. Walbum, B. Brandstrup, A. Celis, M. Puype, J. Van Damme and J. Vandekerckhove A comprehensive two-dimensional gel protein database of noncultured unfractio-J. E. Celis, P. Madsen, nated normal human epidermal keratinocytes: Towards an integrated approach to H. H. Rasmussen, H. Leffers, the study of cell proliferation, differentiation and skin diseases B. Honoré, B. Gesser, K. Dejgaard, E. Olsen, N. Magnusson, J. Kiil, A. Celis, J. B. Lauridsen, B. Basse, G. P. Ratz, A. H. Andersen, E. Walbum, B. Brandstrup, P. S. Pedersen, N. J. Brandt, M. Puype, J. Van Damme and J. Vandekerckhove Microsequencing of proteins recorded in human two-dimensional gel protein H. H. Rasmussen, J. Van Damme, 873 A two-dimensional gel database of human plasma proteins

A two-dimensional gel database of rat liver proteins useful in gene regulation and drug effects studies M. Puype, B. Gesser, J. E. Celis and J. Vandekerckhove 883 N. L. Anderson and N. G. Anderson N. L. Anderson, R. Esquer-Blasco, 907 J.-P. Hofmann and N. G. Anderson The rat liver epithelial (RLE) cell protein database 931 P. J. Wirth, L.-di Luo, Y. Fujimoto, H. C. Bisgaard and A. D. Olson The gene-protein database of Escherichia coli: Edition 4 955 R. A. VanBogelen and F. C. Neidhardt Miscellaneous 995

For submission of papers, see Instructions to Authors (last page of this issue)